

NOV 5 2002

K023197

510(K) SUMMARY

Submitter of 510(k): SB LUCIUS, INC.
9778 Katella Ave. Ste.205, Anaheim, CA 92804
Phone: (714) 530-2814, Fax : (714) 530-3448,

Contact person: Dae Kyu Chang
Phone: (714) 530-2814
Fax: (714) 530-3448
E-mail: Gatorkr@yahoo.co.kr

Date of Summary: Sept. 16, 2002

Trade name: CERANUM 82
Common: Dental casting alloy
Classification name: Gold based alloys and precious metal alloys for clinical use
Product code: EJT
Classification: Class II

Legally marketed device: Metalor's V-Supragold
510(k) number: K993508

SUMMARY OF THE TECHNOLOGICAL CHARACTERISTICS

Test methods applied: as in ANSI/ADA 5 and ISO 9693

Comparison of composition:

COMPOSITION (WEIGHT, %)

Device Name	Au (%)	Pt (%)	Pd (%)	ln (%)	Ir	Ta
V-Supragold	83.00	8.8	5.20	2.8	<1.0%	<1.0%
Ceranium 82	82.00	8.87	3.8	2.0	0.1	0.25

Comparison of physical and mechanical properties:

Alloy	Melting Point Range (°F)	Hardness (Vickers)	Yield Strength (MPa)	Elongation (%)	CTE (x10 ⁻⁶ /°C)	Density (g/cm ³)
V-Supragold	2,012-2,210	190	425	7.0	14.3	
Ceranium 82	2,192-2,372	165	320	9.3	14.1	19.2

Discussion:

Since the composition of the legally marketed alloy and the new device is very similar, it may be assumed that the biological compatibility of the alloys is also very similar.

Conclusion:

The main elements and their concentration are almost identical. CERANUM 82 is a High Noble, Micro-fine, Yellow, Gold based alloy to be used for single crowns, bridge, implant superstructures and substrate for medium expansion with lower fusing porcelains and indirect restorative composites. CERANUM 82 is a high gold ceramic alloy, which heightens the porcelain esthetics of the restoration and provides strength, durability and color of gold. Despite minor differences in the materials, we believe that CERANUM 82 is a substantially equivalent to Metalor's V-Supragold. These changes do not affect safety or effectiveness.



Food and Drug Administration
9200 Corporate Boulevard
Rockville MD 20850

NOV 5 2002

Mr. Dae-Kyu Chang
SB Lucius, Incorporated
9778 Katella Avenue, Suite 205
Anaheim, California 92804

Re: K023197
Trade/Device Name: Ceranum 82
Regulation Number: 21 CFR 872.3060
Regulation Name: Gold-Based Alloys and Precious Metal Alloys for Clinical Use
Regulatory Class: II
Product Code: EJT
Dated: September 16, 2002
Received: September 25, 2002

Dear Mr. Chang:

We have reviewed your Section 510(k) premarket notification of intent to market the device referenced above and have determined the device is substantially equivalent (for the indications for use stated in the enclosure) to legally marketed predicate devices marketed in interstate commerce prior to May 28, 1976, the enactment date of the Medical Device Amendments, or to devices that have been reclassified in accordance with the provisions of the Federal Food, Drug, and Cosmetic Act (Act) that do not require approval of a premarket approval application (PMA). You may, therefore, market the device, subject to the general controls provisions of the Act. The general controls provisions of the Act include requirements for annual registration, listing of devices, good manufacturing practice, labeling, and prohibitions against misbranding and adulteration.

If your device is classified (see above) into either class II (Special Controls) or class III (PMA), it may be subject to such additional controls. Existing major regulations affecting your device can be found in the Code of Federal Regulations, Title 21, Parts 800 to 898. In addition, FDA may publish further announcements concerning your device in the Federal Register.

Please be advised that FDA's issuance of a substantial equivalence determination does not mean that FDA has made a determination that your device complies with other requirements of the Act or any Federal statutes and regulations administered by other Federal agencies.


Page 2 -- Mr. Chang

You must comply with all the Act's requirements, including, but not limited to: registration and listing (21 CFR Part 807); labeling (21 CFR Part 801); good manufacturing practice requirements as set forth in the quality systems (QS) regulation (21 CFR Part 820); and if applicable, the electronic product radiation control provisions (Sections 531-542 of the Act); 21 CFR 1000-1050.

This letter will allow you to begin marketing your device as described in your Section 510(k) premarket notification. The FDA finding of substantial equivalence of your device to a legally marketed predicate device results in a classification for your device and thus, permits your device to proceed to the market.

If you desire specific advice for your device on our labeling regulation (21 CFR Part 801 and additionally 21 CFR Part 809.10 for in vitro diagnostic devices), please contact the Office of Compliance at (301) 594-4613. Additionally, for questions on the promotion and advertising of your device, please contact the Office of Compliance at (301) 594-4639. Also, please note the regulation entitled, "Misbranding by reference to premarket notification" (21CFR Part 807.97). Other general information on your responsibilities under the Act may be obtained from the Division of Small Manufacturers, International and Consumer Assistance at its toll-free number (800) 638-2041 or (301) 443-6597 or at its Internet address <http://www.fda.gov/cdrh/dsma/dsmamain.html>

Sincerely yours,

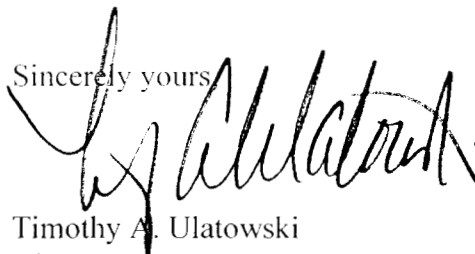

Timothy A. Ulatowski
Director

Division of Anesthesiology, General Hospital,
Infection Control and Dental Devices
Office of Device Evaluation
Center for Devices and
Radiological Health

Enclosure

SB LUCIUS, INC.

9778 Katella Ave. Ste. 205, Anaheim, CA 92804

Phone: (714) 530-2814,
Fax: (714) 530-3448.

INDICATIONS FOR USE

510(K) Number : K023197

Device Name(s) : CERANUM 82

CERANUM 82 is intended for manufacturing:

- Single Crown
- Bridges
- Implant superstructures

(PLEASE DO NOT WRITE BELOW THIS LINE CONTINUE ON ANOTHER
PAGE IF NEEDED)

CONCURRENCE OF CHRD, OFFICE OF DEVICE EVALUATION(OED)

AS Per DMS for Dr. S. Runner

(Division Sign-Off)
Division of Anesthesiology, General Hospital,
Infection Control, Dental Devices

510(k) Number K023197