

MAY 16 2003

K023908  
page 1 of 1

**510(k) SUMMARY**

**SUBMITTED BY**

Prosie Rey-Fessler, RAC  
Director, Regulatory Affairs and Quality Assurance  
INTERPORE CROSS International  
181 Technology Drive  
Irvine, California 92618

**CLASSIFICATION, COMMON OR USUAL NAME, DEVICE NAME**

Classification Name: Cement Restrictor  
Common/Usual Name: Cement Restrictor  
Product Classification: Unclassified  
Proprietary Name: Interpore Cross CEMENT RESTRICTOR

**PREDICATE DEVICES**

Medtronic Sofamor Danek Cement Restrictor  
Macropore *IB* Resorbable Plug  
Pro Osteon 500R Resorbable Bone Graft Substitute  
Pro Osteon 500R Resorbable Bone Void Filler

**INDICATIONS-FOR-USE**

The Interpore Cross CEMENT RESTRICTOR is intended for use as a cement restrictor in the femur, tibia and/or humerus.

**DEVICE DESCRIPTION**

The Interpore Cross CEMENT RESTRICTOR is a resorbable implant manufactured from a polymer-ceramic composite. The Interpore Cross CEMENT RESTRICTOR is designed to wedge into the medullary canal during joint arthroplasty to prevent flow of cement into the canal and allow for pressurized filling of the implantation site. It is provided in various shapes and sizes as needed for particular surgical procedures.

**COMPARISON TO THE PREDICATE DEVICE**

The Interpore Cross CEMENT RESTRICTOR is substantially equivalent to the cited predicate devices based on the indications for use, design features, principles of operation and materials.



DEPARTMENT OF HEALTH & HUMAN SERVICES

Public Health Service

Food and Drug Administration  
9200 Corporate Boulevard  
Rockville MD 20850

Ms. Prosie Rey-Fessler, RAC  
Director, Regulatory Affairs and Quality Assurance  
INTERPORE CROSS International  
181 Technology Drive  
Irvine, California 92618

MAY 16 2003

Re: K023908  
Trade/Device Name: Interpore Cross Cement Restrictor  
Regulation Number: 21 CFR 878.3300  
Regulation Name: Surgical mesh  
Regulatory Class: II  
Product Code: JDK  
Dated: February 28, 2003  
Received: March 3, 2003

Dear Ms. Rey-Fessler:

We have reviewed your Section 510(k) premarket notification of intent to market the device referenced above and we have determined the device is substantially equivalent (for the indications for use stated in the enclosure) to legally marketed predicate devices marketed in interstate commerce prior to May 28, 1976, the enactment date of the Medical Device Amendments, or to devices that have been reclassified in accordance with the provisions of the Federal Food, Drug, and Cosmetic Act (Act) that do not require approval of a premarket approval (PMA) application. You may, therefore, market the device, subject to the general controls provisions of the Act and the limitations described below. The general controls provisions of the Act include requirements for annual registration, listing of devices, good manufacturing practice, labeling, and prohibitions against misbranding and adulteration.

The Office of Device Evaluation has determined that there is a reasonable likelihood that this device will be used for an intended use not identified in the proposed labeling and that such use could cause harm. Therefore, in accordance with Section 513(i)(1)(E) of the Act, the following limitation must appear in the Warnings section of the device's package insert and also as a Warning on the product label:

WARNING: THIS DEVICE IS NOT INTENDED FOR ANY SPINAL INDICATIONS.

THE SAFETY AND EFFECTIVENESS OF THIS DEVICE WHEN  
IMPLANTED IN THE SPINE HAVE NOT BEEN ESTABLISHED.

If your device is classified (see above) into either class II (Special Controls) or class III (PMA), it may be subject to such additional controls. Existing major regulations affecting your device can

be found in the Code of Federal Regulations, Title 21, Parts 800 to 898. In addition, FDA may publish further announcements concerning your device in the Federal Register.

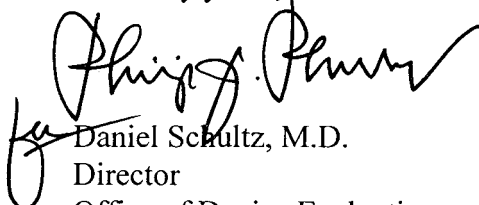
Please be advised that FDA's issuance of a substantial equivalence determination does not mean that FDA has made a determination that your device complies with other requirements of the Act or any Federal statutes and regulations administered by other Federal agencies. You must comply with all the Act's requirements, including, but not limited to: registration and listing (21 CFR Part 807); labeling (21 CFR Part 801); good manufacturing practice requirements as set forth in the quality systems (QS) regulation (21 CFR Part 820); and if applicable, the electronic product radiation control provisions (Sections 531-542 of the Act); 21 CFR 1000-1050.

The FDA finding of substantial equivalence of your device to a legally marketed predicate device results in a classification for your device and permits your device to proceed to the market. This letter will allow you to begin marketing your device as described in your 510(k) premarket notification if the limitation statement above is added to your labeling, as described.

Please note that the above labeling limitations are required by Section 513(i)(1)(E) of the Act. Therefore, a new 510(k) is required before these limitations are modified in any way or removed from the device's labeling.

If you desire specific advice for your device on our labeling regulation (21 CFR Part 801 and additionally 21 CFR Part 809.10 for in vitro diagnostic devices), please contact the Office of Compliance at (301) 594-4659. Additionally, for questions on the promotion and advertising of your device, please contact the Office of Compliance at (301) 594-4639. Also, please note the regulation entitled, "Misbranding by reference to premarket notification" (21 CFR Part 807.97). Other general information on your responsibilities under the Act may be obtained from the Division of Small Manufacturers, International and Consumer Assistance at its toll-free number (800) 638-2041 or (301) 443-6597 or at its Internet address: <http://www.fda.gov/cdrh/dsma/dsmamain.html>

Sincerely yours,

  
Daniel Schultz, M.D.  
Director  
Office of Device Evaluation  
Center for Devices and  
Radiological Health

Enclosure

510(k) Number (if known): K023908

Device Name: INTERPORE CROSS CEMENT RESTRICTOR

**Indications-For-Use:**

The Interpore Cross CEMENT RESTRICTOR is indicated for use as a cement restrictor in the femur, tibia and/or humerus.

(PLEASE DO NOT WRITE BELOW THIS LINE – CONTINUE ON ANOTHER PAGE IF NEEDED)

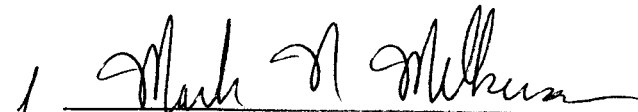
Concurrence of CDRH, Office of Device Evaluation (ODE)

Prescription Use \_\_\_\_\_  
(PER 21 CFR 801.109)

OR

Over-The-Counter Use

(Optional Format 1-2-96)

  
(Division Signature)  
Division of General Restorative  
and Neurological Devices

510(k) Number K023908