

510(k) Summary of Safety and Effectiveness

K030893

LEVO *combi*

APR - 2 2003

Submitter: LEVO AG
P.O. Box 1272
Anglikerstrasse 20
CH-5610 Wohlen, Switzerland

Phone: +41 56-618 44 11
Fax: +41 56-618 44 10

Contact Person: Mr. Thomas Raeber, Managing Director or
Dr. Kurt H. Fischer, CEO

E-Mail-Address: raeber@levo.ch
fischer@levo.ch

Device Name: LEVO *combi*

Classification Name: Electrical Power Wheelchair, Stand-up

Identification of Predicate Devices: LEVO *mobil* LCM (K 963817)

Intended Use: The LEVO *combi* stand-up electrical power wheelchair is a product which changes peoples position not only from sitting to standing and standing to sitting but also reclines and lifts the seat and back position. The product provides indoor and outdoor mobility.

Description of the Device: The LEVO *combi* powered wheelchair is centre wheel driven, battery powered, motor driven and is controlled by the Penny & Giles compact power wheelchair controller Pilot+. The joystick is integrated in the controller.
The wheelchair is powered by two 12V55Ah batteries with a theoretical driving range of 35km on the fully charged batteries.
The wheelchair consists of the following basic sub-sections:

- Base with two direct-drive units with integrated parking brakes, two 12V55Ah batteries, two 14"x2¾" centre wheels, two 8"x2" front wheels, and one 5½"x1½" rear wheel with suspension.
- Penny & Giles power electronics Pilot+
- Seating system including mechanisms for all possible positions, stand-up actuator and optional lift and/or backrest actuator.

The base is of welded steel construction and includes the base frame, front castor wheels, centre driving wheels with drive unit (motor/gear/brakes), batteries and pivoting rear castor. The motor controller is mounted to the left or right armrest, depending on the user's needs.

Substantial
Equivalence:

The **LEVO combi** is substantial equivalent with respect to intended use, and materials to the predicate device.

The key changes are rear suspension, outdoor mobility, the use of the Pilot+ controller from Penny & Giles instead of the PG8 controller from Penny & Giles, as well as additional possibilities to change positions (reclining up to lying, relax position).

Safety and
Effectiveness:

The **LEVO combi** was basically developed on experience of the **LEVO mobil** LCM. However being able to provide additional outdoor mobility as well as further improvements in the geometric for an excellent biomechanical response to the user's body, there are a few additional functions and adjustments integrated. The rest of it is simply colour and design changes. So the **LEVO combi** has in substantial the same technological characteristics and the same safety and effectiveness as the predicate device(s) and the minor changes declared in the submission do not raise new questions of safety and effectiveness.



DEPARTMENT OF HEALTH & HUMAN SERVICES

Public Health Service

Food and Drug Administration
9200 Corporate Boulevard
Rockville MD 20850

APR - 2 2003

LEVO AG
c/o Mr. Stefan Preiss
TUV America, Inc.
1775 Old Highway 8
New Brighton, MN 55112

Re: K030893
Trade/Device Name: LEVO combi
Regulation Number: 21 CFR 890.3900
Regulation Name: Standup wheelchair
Regulatory Class: II
Product Code: IPL
Dated: March 21, 2003
Received: March 21, 2003

Dear Mr. Preiss:

We have reviewed your Section 510(k) premarket notification of intent to market the device referenced above and have determined the device is substantially equivalent (for the indications for use stated in the enclosure) to legally marketed predicate devices marketed in interstate commerce prior to May 28, 1976, the enactment date of the Medical Device Amendments, or to devices that have been reclassified in accordance with the provisions of the Federal Food, Drug, and Cosmetic Act (Act) that do not require approval of a premarket approval application (PMA). You may, therefore, market the device, subject to the general controls provisions of the Act. The general controls provisions of the Act include requirements for annual registration, listing of devices, good manufacturing practice, labeling, and prohibitions against misbranding and adulteration.


If your device is classified (see above) into either class II (Special Controls) or class III (PMA), it may be subject to such additional controls. Existing major regulations affecting your device can be found in the Code of Federal Regulations, Title 21, Parts 800 to 898. In addition, FDA may publish further announcements concerning your device in the Federal Register.

Please be advised that FDA's issuance of a substantial equivalence determination does not mean that FDA has made a determination that your device complies with other requirements of the Act or any Federal statutes and regulations administered by other Federal agencies. You must comply with all the Act's requirements, including, but not limited to: registration and listing (21 CFR Part 807); labeling (21 CFR Part 801); good manufacturing practice requirements as set forth in the quality systems (QS) regulation (21 CFR Part 820); and if applicable, the electronic product radiation control provisions (Sections 531-542 of the Act); 21 CFR 1000-1050.

This letter will allow you to begin marketing your device as described in your Section 510(k) premarket notification. The FDA finding of substantial equivalence of your device to a legally marketed predicate device results in a classification for your device and thus, permits your device to proceed to the market.

If you desire specific advice for your device on our labeling regulation (21 CFR Part 801), please contact the Office of Compliance at (301) 594-4659. Also, please note the regulation entitled, "Misbranding by reference to premarket notification" (21CFR Part 807.97). You may obtain other general information on your responsibilities under the Act from the Division of Small Manufacturers, International and Consumer Assistance at its toll-free number (800) 638-2041 or (301) 443-6597 or at its Internet address <http://www.fda.gov/cdrh/dsma/dsmamain.html>

Sincerely yours,

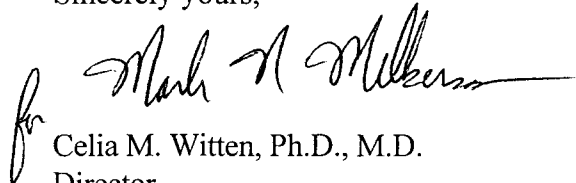for Celia M. Witten, Ph.D., M.D.

Director
Division of General, Restorative
and Neurological Devices
Office of Device Evaluation
Center for Devices and
Radiological Health

Enclosure

Statement of Indication for Use


Device Name: **LEVO combi**

Intended Use:

The **LEVO combi** stand-up electrical power wheelchair is a product which changes peoples position not only from sitting to standing and standing to sitting but also reclines and lifts the seat and back position. The product provides indoor and outdoor mobility.

Target Population:

For all individuals who need a power wheelchair with the possibility to change positions and who can not stand on their feet themselves such as people with spinal cord injury, spina bifida, cerebral palsy, multiple sclerosis, muscular dystrophy, polio, rheumatism, etc..

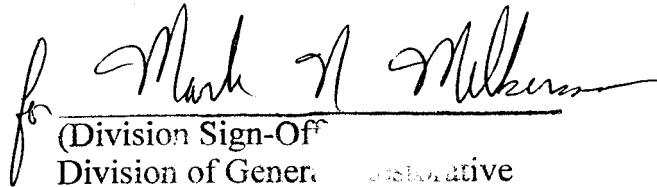

(Division Sign-Off)

Division of General Rehabilitative
and Neurological Services

510(k) Number _____

K030893

PM
T