



Food and Drug Administration
9200 Corporate Boulevard
Rockville MD 20850

MAY 4 - 2004

Asia Pacific Microsystems, Inc.
c/o Ke-Min Jen, Ph.D.
ROC Chinese-European Industrial Research Society
No 2, R&D Road 6, Hsin Chu
Science-Based Industrial Park
Hsin Chu, (TAIWAN) 300

Re: K040159
APM Blood Pressure Monitor, BP108A
Regulation Number: 21 CFR 870.1130
Regulation Name: Noninvasive Blood Pressure Measurement System
Regulatory Class: Class II (two)
Product Code: DXN
Dated: January 10, 2004
Received: January 23, 2004

Dear Dr. Jen:

This letter corrects our substantially equivalent letter of March 22, 2004 regarding the incorrect spelling of the name of your company.

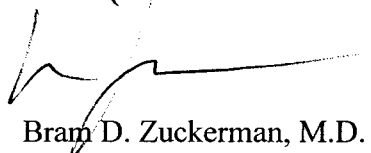
We have reviewed your Section 510(k) premarket notification of intent to market the device referenced above and have determined the device is substantially equivalent (for the indications for use stated in the enclosure) to legally marketed predicate devices marketed in interstate commerce prior to May 28, 1976, the enactment date of the Medical Device Amendments, or to devices that have been reclassified in accordance with the provisions of the Federal Food, Drug, and Cosmetic Act (Act) that do not require approval of a premarket approval application (PMA). You may, therefore, market the device, subject to the general controls provisions of the Act. The general controls provisions of the Act include requirements for annual registration, listing of devices, good manufacturing practice, labeling, and prohibitions against misbranding and adulteration.

If your device is classified (see above) into either class II (Special Controls) or class III (PMA), it may be subject to such additional controls. Existing major regulations affecting your device can be found in the Code of Federal Regulations, Title 21, Parts 800 to 898. In addition, FDA may publish further announcements concerning your device in the Federal Register.

Please be advised that FDA's issuance of a substantial equivalence determination does not mean that FDA has made a determination that your device complies with other requirements of the Act or any Federal statutes and regulations administered by other Federal agencies. You must comply with all the Act's requirements, including, but not limited to: registration and listing (21 CFR Part 807); labeling (21 CFR Part 801); good manufacturing practice requirements as set forth in the quality systems (QS) regulation (21 CFR Part 820); and if applicable, the electronic product radiation control provisions (Sections 531-542 of the Act); 21 CFR 1000-1050. This letter will allow you to begin marketing your device as described in your Section 510(k) premarket notification. The FDA finding of substantial equivalence of your device to a legally marketed predicate device results in a classification for your device and thus, permits your device to proceed to the market.

If you desire specific advice for your device on our labeling regulation (21 CFR Part 801), please contact the Office of Compliance at (301) 594-4648. Also, please note the regulation entitled, "Misbranding by reference to premarket notification" (21CFR Part 807.97). You may obtain other general information on your responsibilities under the Act from the Division of Small Manufacturers, International and Consumer Assistance at its toll-free number (800) 638-2041 or (301) 443-6597 or at its Internet address <http://www.fda.gov/cdrh/dsma/dsmamain.html>

Sincerely yours,

A handwritten signature in black ink, appearing to read 'Bram D. Zuckerman', is written over a horizontal line.

Bram D. Zuckerman, M.D.
Director
Division of Cardiovascular Devices
Office of Device Evaluation
Center for Devices and
Radiological Health

Enclosure

Indications for Use

510(k) Number (if known): K040159

Device Name: APM Blood Pressure Monitor, BP108A

Indications For Use:

The device as a noninvasive blood pressure measurement system intended to measure the systolic and diastolic blood pressures and pulse rate on an adult individual, over age 18, at home by using a non-invasive technique in which an inflatable cuff is wrapped around the wrist. The cuff circumference is limited to 5.25"-7.75".

Prescription Use _____
(Part 21 CFR 801 Subpart D)

AND/OR

Over-The-Counter Use ✓
(21 CFR 807 Subpart C)

(PLEASE DO NOT WRITE BELOW THIS LINE-CONTINUE ON ANOTHER PAGE IF NEEDED)

Concurrence of CDRH, Office of Device Evaluation (ODE)

Dan R. Lockner
(Division Sign-Off)
Division of Cardiovascular Devices

510(k) Number K040159

Page 1 of 1

3

K040159

p1/e-

Asia Pacific Microsystems, Inc

No.2,R&D Road 6, Hsinchu Science-based Industrial Park

Hsin Chu, Taiwan, R.O.C

MAR 22 2004

Tel : +886-3-666-1188 Fax : +886-3-666-1199

e-mail : mslin@apmsinc.com

http:// www.apmsinc.com

“ 510(k) SUMMARY ”

Submitter's Name: **Asia Pacific Microsystems, Inc.**

No.2,R&D Road 6, Hsinchu Science-based Industrial Park

Hsin Chu, Taiwan, R.O.C

Telephone: 886-3-666-1188 Fax: 886-3-666-1199

e-mail : mslin@apmsinc.com

Date summary prepared:

January 10, 2004

Device Name:

Proprietary Name: APM Blood Pressure Monitor, BP108A

Common or Usual Name: NONINVASIVE BLOOD PRESSURE MEASUREMENT SYSTEM

Classification Name: Blood Pressure Monitor, Class II,

~~21~~ CFR 870.1130

Indications for Use:

The device is a noninvasive blood pressure measurement system intended to measure the systolic and diastolic blood pressures and pulse rate of an adult individual, over age 18, at home by using a non-invasive technique in which an inflatable cuff is wrapped around the wrist. The cuff circumference is limited to be 5.25" – 7.75".

Description of the device:

APM BP108A uses the Oscillometric method to measure the blood pressure. The Oscillometric method is adopted clinically to measure the blood pressure recently. It is not needed to use the stethoscope, as in the traditional measuring method, to monitor the Korotkov sound when deciding the systolic or diastolic pressure. The Oscillometric method senses the vibrating signal via the closed air pipe system and utilizes the microcomputer to automatically sense the characteristics of the pulse signal. Through simple calculation, clinically proven, the reading can reflect the accurate real blood pressure, and the systolic pressure is defined as the pressure when the cuff pressure oscillating amplitude begins to increase and the diastolic pressure as the pressure when the cuff pressure oscillating amplitude stops decreasing.

Performance Testing:

Electric Safety Requirement Test Report of EN 60601-1:1990 & EN 1060-1/
EN1060-3, and EMC test report of EN 60601-1-2 (EN 55011:1991 and EN
61000-4-2:1995)

ANSI/AAMI SP10-1992 Electronic or Automated Sphygmomanometers

Legally marketed device for substantial equivalence comparison:

Eikon Automatic Digital Blood Pressure Monitor, HD-400M (**K021239**)

Summary for substantial equivalence comparison:

Same characteristics: intended use, technological characteristics, power supply, display,
measuring range, accuracy, operating and storage environments.

Different characteristics: memory, dimensions, weight, and storage environments.

1. As we can understand, the memory feature is to memorize the measurement data taken previously and is related to the usage convenience, not to raise any safety or effectiveness hazard.
2. The differences between dimensions and weight are related to the designing aspects.
These differences are not to raise any safety or effectiveness aspect.

They are decided to be substantially equivalent.