

K042009

Special 510(k) Premarket Notification
USA Instruments – Mark III Shoulder Coil
July 23, 2004



USA Instruments, Inc.

AUG 11 2004

SUMMARY OF SAFETY AND EFFECTIVENESS

1. Device Name : Magnetic Resonance Imaging Accessory
2. Proprietary Name : Mark III Phased Array Shoulder Coil
3. Classification : Class II
4. Establishment Registration #: 1529041
5. Manufacture Facility Location: USA Instruments, Inc., 1515 Danner Drive,
Aurora, Ohio 44202, USA
Telephone: 330-995-8500; Fax: 330-562-1422.
6. Performance Standard: No applicable performance standards have been issued
under Section 514 of the Food, Drug and Cosmetic Act.
7. Intended Use: The Mark III Phased Array Shoulder Coil is a receive-only phased
array RF coil, used for obtaining diagnostic images of the shoulder
and adjacent regions in Magnetic Resonance Imaging systems. The
indications for use are the same as for standard MR Imaging. The
Mark III Phased Array Coil is designed for use with the Signa
EXCITE 3T MRI scanner manufactured by GE Medical Systems,
Inc.
8. Device Description : The Mark III Phased Array Coil consists of three volume RF coil
elements in a phased array design. The coil elements and
associated circuitry are enclosed to prevent any exposure to the
patient or environment. The coil electronics are enclosed in both
the rigid housing and the vinyl coated PVC foam. The coil is
positioned on the patient's shoulder during imaging.

9. Safety and Effectiveness:

Mark III Phased Array Shoulder Coil Product Features	Comparison to Predicate or other 510(k) cleared products
Intended Use Shoulder imaging applications	-Similar to the Mark 9000 Shoulder Coil manufactured by USA Instruments Inc. (K010946) -Similar to the Phased Array Shoulder Coil manufactured by Medical Advances Inc. (K945778)
Indications for Use Identical to routine MRI imaging	-Similar to the Mark 9000 Shoulder Coil manufactured by USA Instruments Inc. (K010946) -Similar to the Phased Array Shoulder Coil manufactured by Medical Advances Inc. (K945778)
Coil Enclosure Material Vinyl coated PVC foam Flame retardent Polyurethane Plastic PVC Plastic	-Similar to the Mark 9000 Shoulder Coil manufactured by USA Instruments Inc. (K010946) -Similar to the Magna 5000 Phased Array CTL Spine Coil manufactured by USA Instruments, Inc. (K994345)
Coil Design Receive-only phased array design	-Similar to the Mark 9000 Shoulder Coil manufactured by USA Instruments Inc. (K010946)
Decoupling Switching diode decoupling	-Similar to the Mark 9000 Shoulder Coil manufactured by USA Instruments Inc. (K010946)
Prevention of RF Burns Does not transmit RF power; decoupling isolates the coil elements from RF fields during RF transmission; coil elements and circuitry are enclosed in a non-conductive housing.	-Similar to the Mark 9000 Shoulder Coil manufactured by USA Instruments Inc. (K010946) -Similar to the Excalibur 9000 TotalSENSE™ Quad Array Head Coil manufactured by USA Instruments, Inc. (K022582)
Radio Frequency Absorption Coil is a receive only coil and does not transmit RF power; power deposition during imaging is limited by SAR algorithm	-Similar to the Mark 9000 Shoulder Coil manufactured by USA Instruments Inc. (K010946)
Formation of Resonant Loop Decoupling isolates the coil elements from RF fields during RF transmission; length of cable and stiffness does not permit looping	-Similar to the Mark 9000 Shoulder Coil manufactured by USA Instruments Inc. (K010946)



Food and Drug Administration
9200 Corporate Boulevard
Rockville MD 20850

AUG 11 2004

Mr. James Wrenn
QA and Regulatory Manager
USA Instruments, Inc.
1515 Danner Drive
AURORA OH 44202

Re: K042009
Trade/Device Name: Mark III Phased Array
Shoulder Coil
Regulation Number: 21 CFR 892.1000
Regulation Name: Magnetic resonance
diagnostic device
Regulatory Class: II
Product Code: 90 MOS
Dated: July 23, 2004
Received: July 26, 2004

Dear Mr. Wrenn:

We have reviewed your Section 510(k) premarket notification of intent to market the device referenced above and have determined the device is substantially equivalent (for the indications for use stated in the enclosure) to legally marketed predicate devices marketed in interstate commerce prior to May 28, 1976, the enactment date of the Medical Device Amendments, or to devices that have been reclassified in accordance with the provisions of the Federal Food, Drug, and Cosmetic Act (Act) that do not require approval of a premarket approval application (PMA). You may, therefore, market the device, subject to the general controls provisions of the Act. The general controls provisions of the Act include requirements for annual registration, listing of devices, good manufacturing practice, labeling, and prohibitions against misbranding and adulteration.

If your device is classified (see above) into either class II (Special Controls) or class III (PMA), it may be subject to such additional controls. Existing major regulations affecting your device can be found in the Code of Federal Regulations, Title 21, Parts 800 to 898. In addition, FDA may publish further announcements concerning your device in the Federal Register.

Please be advised that FDA's issuance of a substantial equivalence determination does not mean that FDA has made a determination that your device complies with other requirements of the Act or any Federal statutes and regulations administered by other Federal agencies. You must comply with all the Act's requirements, including, but not limited to: registration and listing (21 CFR Part 807); labeling (21 CFR Part 801); good manufacturing practice requirements as set forth in the quality systems (QS) regulation (21 CFR Part 820); and if applicable, the electronic product radiation control provisions (Sections 531-542 of the Act); 21 CFR 1000-1050.


This letter will allow you to begin marketing your device as described in your Section 510(k) premarket notification. The FDA finding of substantial equivalence of your device to a legally marketed predicate device results in a classification for your device and thus, permits your device to proceed to the market.

If you desire specific advice for your device on our labeling regulation (21 CFR Part 801), please contact the Office of Compliance at one of the following numbers, based on the regulation number at the top of the letter:

8xx.1xxx	(301) 594-4591
876.2xxx, 3xxx, 4xxx, 5xxx	(301) 594-4616
884.2xxx, 3xxx, 4xxx, 5xxx, 6xxx	(301) 594-4616
892.2xxx, 3xxx, 4xxx, 5xxx	(301) 594-4654
Other	(301) 594-4692

Additionally, for questions on the promotion and advertising of your device, please contact the Office of Compliance at (301) 594-4639. Also, please note the regulation entitled, "Misbranding by reference to premarket notification" (21CFR Part 807.97) you may obtain. Other general information on your responsibilities under the Act may be obtained from the Division of Small Manufacturers, International and Consumer Assistance at its toll-free number (800) 638-2041 or (301) 443-6597 or at its Internet address <http://www.fda.gov/cdrh/dsma/dsmamain.html>.

Sincerely yours,

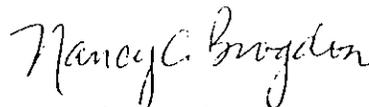

Nancy C. Brogdon
Director, Division of Reproductive,
Abdominal and Radiological Devices
Office of Device Evaluation
Center for Devices and Radiological Health

Enclosure

STATEMENT OF INTENDED USE

510(k) Number (if known): 1K042009

Device Name: Mark III Phased Array Shoulder Coil

Indications for Use:

The Mark III Phased Array Shoulder Coil is a receive-only phased array RF coil, used for obtaining diagnostic images of the shoulder and adjacent regions in Magnetic Resonance Imaging systems. The indications for use are the same as for standard MR Imaging. The Mark III Phased Array Coil is designed for use with the Signa EXCITE 3.0T MRI scanner manufactured by GE Medical Systems, Inc

(PLEASE DO NOT WRITE BELOW THIS LINE – CONTINUE ON ANOTHER PAGE IF NEEDED)

Concurrence of CDRH, Office of Device Evaluation (ODE)

Prescription Use OR Over-The-Counter Use
(Per 21 CFR 801-109)

Nancy Brogdon
(Division Sign-Off)
Division of Reproductive, Abdominal,
and Radiological Devices
510(k) Number K042009