

MAR 3 - 2005

Attachment (D) 510(k) Summary

1. DATE PREPARED

February 4, 2005

2. SPONSOR INFORMATION

A&D Engineering, Inc.
Mr. Jerry Wang
1555 McCandless Drive, Milpitas, CA 95035
Tel: 408-518-5113
Fax: 408-635-2313
Email: jwang@andmedical.com

3. DEVICE NAME

Proprietary Name: A&D Medical UA-767PBT Digital Blood Pressure Monitor

Common/Usual Name: Blood Pressure Monitor

Classification name: Non-invasive blood pressure measurement System
21 CFR 870-1130, Class II, 74DXN.

4. DEVICE DESCRIPTION AND INTENDED USE

The A&D Medical UA-767PBT digital blood pressure monitor is intended for use by adults for measuring the systolic and diastolic blood pressure and pulse rate.

5. PREDICATE DEVCIE

It is substantially equivalent to the following device:

A&D UA-767PC, FDA 510(k) K982481. Issued on Jan. 13, 1999
A&D UA-767BT, FDA 510(k) K040371. Issued on May 19, 2004
A&D UA-787, FDA 510(k) K012472. Issued on Feb. 15, 2002

6. TECHNOLOGICAL CHARACTERISTICS

UA-767PBT uses an inflated cuff which is wrapped around the upper arm. The cuff is inflated by an electrical air pump. The systolic and diastolic blood pressures are determined by oscillometric method. The deflation rate is controlled by a preset mechanical valve at a constant rate. At any moment of measurement, the user can deflate the cuff by pressing the blue "START" button. The measurement results are displayed on the LCD and transmitted to a Bluetooth enabled devices, such as a PC, a PDA, a printer, or an access point. UA-767PBT measures blood pressure and pulse rate even when an irregular heartbeat occurs.

7. DEVICE TESTING

A&D Medical UA-767BT digital blood pressure monitor meets NIST/AAMI SP-10 standard and FDA guidance "Non-invasive Blood Pressure (NIBP) Monitor Guidance". Please refer to the table below for the list of AAMI SP-10 tests. It uses the identical software codes and pressure detection related hardware as the predicate devices to determine systolic, diastolic, and pulse rate.

SP-10 Section #	Section Title	Test Results & Comments
4.1.1	General	Conformed
4.1.2.1	Device labeling	Conformed
4.1.2.2	Outer container	Conformed
4.1.3	Information manual	Conformed
4.1.4.1	Component replacement	Conformed
4.1.4.2	Power system labeling	Conformed
4.1.4.3	Labeling for battery-powered devices	Conformed
4.2.1	Storage conditions	Conformed
4.2.2	Operating conditions	Conformed
4.2.3	Vibration and shock	Conformed
4.2.4.1	Voltage range	Conformed
4.2.4.2	Life	Conformed
4.3.1.1	Maximum cuff pressure	Conformed
4.3.1.2	Cuff deflation	Conformed
4.3.2	Electrical safety	Conformed
4.3.3	Conductive components	Conformed
4.4.1	Pressure indicator accuracy	Conformed
4.4.2	Overall system efficacy	Conformed
4.4.2.1	Auscultatory method as the reference standard	Conformed
4.4.2.2	Intra-arterial method as the reference standard	Not applicable
4.4.3	Battery-powered devices	Conformed
4.5	Requirements for devices with manual inflation systems	Conformed



DEPARTMENT OF HEALTH & HUMAN SERVICES

Public Health Service

MAR 3 - 2005

Food and Drug Administration
9200 Corporate Boulevard
Rockville MD 20850

A&D Engineering, Inc.
c/o Mr Jerry Wang
Director of Engineering & QA
1555 McCandless Drive
Milpitas, CA 95025

Re: K043217

Trade Name: UA-767PBT Digital Blood Pressure Monitor
Regulation Number: 21 CFR 870.1130
Regulation Name: Noninvasive Blood Pressure Measurement System
Regulatory Class: Class II (two)
Product Code: DXN
Dated: February 04, 2005
Received: February 07, 2005

Dear Mr. Wang:

We have reviewed your Section 510(k) premarket notification of intent to market the device referenced above and have determined the device is substantially equivalent (for the indications for use stated in the enclosure) to legally marketed predicate devices marketed in interstate commerce prior to May 28, 1976, the enactment date of the Medical Device Amendments, or to devices that have been reclassified in accordance with the provisions of the Federal Food, Drug, and Cosmetic Act (Act) that do not require approval of a premarket approval application (PMA). You may, therefore, market the device, subject to the general controls provisions of the Act. The general controls provisions of the Act include requirements for annual registration, listing of devices, good manufacturing practice, labeling, and prohibitions against misbranding and adulteration.

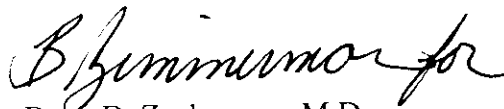
If your device is classified (see above) into either class II (Special Controls) or class III (PMA), it may be subject to such additional controls. Existing major regulations affecting your device can be found in the Code of Federal Regulations, Title 21, Parts 800 to 898. In addition, FDA may publish further announcements concerning your device in the Federal Register.

Page 2 – Mr. Jerry Wang

Please be advised that FDA's issuance of a substantial equivalence determination does not mean that FDA has made a determination that your device complies with other requirements of the Act or any Federal statutes and regulations administered by other Federal agencies. You must comply with all the Act's requirements, including, but not limited to: registration and listing (21 CFR Part 807); labeling (21 CFR Part 801); good manufacturing practice requirements as set forth in the quality systems (QS) regulation (21 CFR Part 820); and if applicable, the electronic product radiation control provisions (Sections 531-542 of the Act); 21 CFR 1000-1050. This letter will allow you to begin marketing your device as described in your Section 510(k) premarket notification. The FDA finding of substantial equivalence of your device to a legally marketed predicate device results in a classification for your device and thus, permits your device to proceed to the market.

If you desire specific advice for your device on our labeling regulation (21 CFR Part 801), please contact the Office of Compliance at (301) 594-4646. Also, please note the regulation entitled, "Misbranding by reference to premarket notification" (21CFR Part 807.97) you may obtain. Other general information on your responsibilities under the Act may be obtained from the Division of Small Manufacturers, International and Consumer Assistance at its toll-free number (800) 638-2041 or (301) 443-6597 or at its Internet address <http://www.fda.gov/cdrh/dsma/dsmamain.html>

Sincerely yours,

A handwritten signature in black ink, appearing to read "B. Zuckerman" with a stylized flourish at the end.

Bram D. Zuckerman, M.D.

Director

Division of Cardiovascular Devices

Office of Device Evaluation

Center for Devices and

Radiological Health

Enclosure

Indications for Use

510(K) Number (if known): K043217

Device Name: A&D Medical UA-767PBT Digital Blood Pressure Monitors

Indications For Use:

The UA-767PBT is designed to be used by end users who are eighteen (18) years and older at home to monitor their blood pressure (systolic and diastolic) and pulse rate. At the end of each measurement, the results will be stored in the UA-767PBT memory. UA-767PBT through its Bluetooth wireless communication port can also transfer the measurements stored in memory to other electronic devices, such as an Access Point, PC, a modem, or a printer. UA-767PBT uses the oscillometric method to conduct the measurement. It is not designed for ambulatory use. The arm circumference range shall be between 5.1 inches (13.0 cm) to 17.7 inches (45.0 cm).

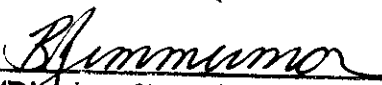
Prescription Use _____
(Part 21 CFR 801 Subpart D)

AND/OR

Over-The-Counter Use X
(21 CFR 807 Subpart C)

(PLEASE DO NOT WRITE BELOW THIS LINE-CONTINUE ON ANOTHER PAGE IF NEEDED)

Concurrence of CDRH, Office of Device Evaluation (ODE)


(Division Sign-Off)
Division of Cardiovascular Devices
510(k) Number K043217