

MAR 25 2005

1650597

Date: March 7<sup>th</sup>, 2005  
Subject: 510(k) Summary of Safety and Effectiveness Information  
for the GE Datex-Ohmeda Engstrom Carestation  
  
Proprietary: GE Datex-Ohmeda Engstrom Carestation  
  
Common: Ventilator, Continuous  
  
Classification: Anesthesiology, 73 CBK, 21 CFR 868.5895

The 510(k) summary of safety and effectiveness information is being submitted in accordance with the requirements of SMDA 1990 and 1992.

The GE Datex-Ohmeda Engstrom Carestation is substantially equivalent to the following currently marketed device:

GE Datex-Ohmeda Engstrom Carestation- Class II - 21CFR868.5895, which has been the subject of a cleared 510(k) with FDA log number K041775.

Siemens Mini C Compressor – Class I – 21CFR 868.6250, which has been the subject of a cleared 510(k) with FDA log number K023354

The Engström Ventilator (EV) is a flexible, adaptable, and intuitive critical care ventilator. A wide selection of performance options gives the user full control of the system configuration. The Engström Carestation is a complete system featuring patient monitoring, patient ventilation, and the capability of interfacing with central information management systems.

The GE Datex-Ohmeda Engstrom Carestation is designed to provide mechanical ventilation for adults and pediatrics weighing 5kg and above having degrees of pulmonary impairment varying from minor to severe. The modes of ventilation are available include:

- Volume Controlled (VCV)
- Pressure Controlled (PCV)
- Pressure Controlled, Volume Guaranteed (PCV-VG)
- Synchronized Intermittent Mandatory Ventilation, Volume Controlled (SIMV-VC)
- Synchronized Intermittent Mandatory Ventilation, Pressure Controlled (SIMV-PC)
- Bi-level Airway Pressure Ventilation
- Constant Positive Airway Pressure/Pressure Support Ventilation (CPAP/PSV)
- Apnea backup (active in Bi-level and CPAP/PSV)

The GE Datex-Ohmeda Engstrom Carestation is a microprocessor based, electronically controlled, pneumatically driven ventilator that includes integrated FiO<sub>2</sub>, airway pressure, spirometry and volume monitoring and an Aerogen Aeroneb Pro nebulzier. Options include integrated respiratory gas monitoring capabilities via various Datex-Ohmeda patient monitoring modules listed in the product labeling and an integrated air compressor.

When supplied as an option, the integrated respiratory gas monitoring is provided via the Datex-Ohmeda M-Gas Module (M-C, M-CO, M-COV, M-COVX, M-CaIO, M-CAiOV, M-CAiOVX,

(rev 3.2 software and higher) K# 001814) or Mini-CO2 Module (K023454) which are physically integrated into the Engstrom Carestation, receive electronic power from the Engstrom Carestation and communicate measured values to the Engstrom Carestation for display on the system display unit.

When supplied as an option, the GE Datex-Ohmeda EV Air Compressor is intended for use as an accessory to provide a dry, filtered, breathable compressed air supply. The compressor has no alarm functions. All alarm functions and reactions to failure of the compressed gas supply, are provided by the Engstrom Carestation as cleared in K041775. The compressor is installed in the base of the ventilator cart. The compressor is powered from AC mains only. A source of compressed oxygen is required to be connected to Engstrom Carestations equipped with the optional compressor.

The system is designed for facility use, including within-facility transport, and should only be used under the orders of a clinician.

The ventilator consists of three main components: a display, a ventilator unit, and an optional module bay. The display allows the user to interface with the system and control settings. The ventilator unit controls electrical power, nebulization, and pneumatic gas flow to and from the patient. The module bay allows the integration of various Datex-Ohmeda patient monitoring modules with the ventilator.

Additional optional accessories include a trolley/cart, compressor, airway modules, module bay, support arm, humidifier and water trap mounting brackets, and auxiliary electrical outlets. The user interface for control of nebulization is provided via the ventilator display unit. The Aerogen Aeroneb Pro Nebulizer (K021175) is provided standard with the unit.

The GE Datex-Ohmeda Engstrom Carestation was designed to comply with the applicable portions of the following voluntary standards;

1. UL 2601 – General requirements for Medical Electrical Equipment
2. ASTM F1100 – Particular Requirements for Critical Care Ventilators
3. EN/IEC 60601-1: General requirements for Medical Electrical Equipment
4. EN/IEC 60601-1-2: 2001 - Medical Electrical Equipment - Electromagnetic Compatibility
5. EN 475 – Electrically Generated Alarm Signals
6. CGA V-1 and ISO 5145 Medical Gas Cylinders – Threaded Cylinders
7. EN 980 Graphical Symbols
8. EN/IEC 60601-2-12, Medical Electrical Equipment – Critical Care Ventilators

The GE Datex-Ohmeda Engstrom Carestation and the currently marketed device are substantially equivalent in design concepts, technologies and materials. The GE Datex-Ohmeda Engstrom Carestation has been validated through rigorous testing that, in part, supports the compliance of GE Datex-Ohmeda Engstrom Carestation to the standards listed above.

Contact: Dan Kosednar, RAC  
Manager, Regulatory Planning and Submissions



DEPARTMENT OF HEALTH & HUMAN SERVICES

Public Health Service

MAR 25 2005

Food and Drug Administration  
9200 Corporate Boulevard  
Rockville MD 20850

Mr. Dan Kosednar  
Manager, Regulatory Planning and Submissions  
Datex-Ohmeda, Incorporated  
Care Business Area  
P.O. Box 7550  
Madison, Wisconsin 53707

Re: K050597

Trade/Device Name: GE Datex-Ohmeda Engstrom Carestation  
Regulation Number: 21 CFR 868.5895  
Regulation Name: Continuous Ventilator  
Regulatory Class: II  
Product Code: CBK  
Dated: March 8, 2005  
Received: March 9, 2005

Dear Mr. Kosednar:

We have reviewed your Section 510(k) premarket notification of intent to market the device referenced above and have determined the device is substantially equivalent (for the indications for use stated in the enclosure) to legally marketed predicate devices marketed in interstate commerce prior to May 28, 1976, the enactment date of the Medical Device Amendments, or to devices that have been reclassified in accordance with the provisions of the Federal Food, Drug, and Cosmetic Act (Act) that do not require approval of a premarket approval application (PMA). You may, therefore, market the device, subject to the general controls provisions of the Act. The general controls provisions of the Act include requirements for annual registration, listing of devices, good manufacturing practice, labeling, and prohibitions against misbranding and adulteration.

If your device is classified (see above) into either class II (Special Controls) or class III (PMA), it may be subject to such additional controls. Existing major regulations affecting your device can be found in the Code of Federal Regulations, Title 21, Parts 800 to 898. In addition, FDA may publish further announcements concerning your device in the Federal Register.

Please be advised that FDA's issuance of a substantial equivalence determination does not mean that FDA has made a determination that your device complies with other requirements of the Act or any Federal statutes and regulations administered by other Federal agencies. You must comply with all the Act's requirements, including, but not limited to: registration and listing (21 CFR Part 807); labeling (21 CFR Part 801); good manufacturing practice requirements as set forth in the quality systems (QS) regulation (21 CFR Part 820); and if applicable, the electronic product radiation control provisions (Sections 531-542 of the Act); 21 CFR 1000-1050.

This letter will allow you to begin marketing your device as described in your Section 510(k) premarket notification. The FDA finding of substantial equivalence of your device to a legally marketed predicate device results in a classification for your device and thus, permits your device to proceed to the market.

If you desire specific advice for your device on our labeling regulation (21 CFR Part 801), please contact the Office of Compliance at (240) 276-0102. Also, please note the regulation entitled, "Misbranding by reference to premarket notification" (21CFR Part 807.97). You may obtain other general information on your responsibilities under the Act from the Division of Small Manufacturers, International and Consumer Assistance at its toll-free number (800) 638-2041 or (301) 443-6597 or at its Internet address <http://www.fda.gov/cdrh/industry/support/index.html>.

Sincerely yours,

A handwritten signature in black ink, appearing to read "Chiu Lin, Ph.D.", written in a cursive style.

Chiu Lin, Ph.D.

Director

Division of Anesthesiology, General Hospital,

Infection Control and Dental Devices

Office of Device Evaluation

Center for Devices and

Radiological Health

Enclosure

## Indications for Use

510(k) Number (if known): K050597

Device Name: GE Datex-Ohmeda Engstrom Carestation

### Indications For Use:

The GE Datex-Ohmeda Engstrom Carestation is designed to provide mechanical ventilation for adults and pediatrics weighing 5kg and above having degrees of pulmonary impairment varying from minor to severe. The modes of ventilation are available include:

- Volume Controlled (VCV)
- Pressure Controlled (PCV)
- Pressure Controlled, Volume Guaranteed (PCV-VG)
- Synchronized Intermittent Mandatory Ventilation, Volume Controlled (SIMV-VC)
- Synchronized Intermittent Mandatory Ventilation, Pressure Controlled (SIMV-PC)
- Bi-level Airway Pressure Ventilation
- Constant Positive Airway Pressure/Pressure Support Ventilation (CPAP/PSV)
- Apnea backup (active in Bi-level and CPAP/PSV)

The GE Datex-Ohmeda Engstrom Carestation is a microprocessor based, electronically controlled, pneumatically driven ventilator that includes integrated FIO<sub>2</sub>, airway pressure, spirometry and volume monitoring and an Aerogen Aeroneb Pro nebulzier. Options include integrated respiratory gas monitoring capabilities via various Datex-Ohmeda patient monitoring modules listed in the product labeling and an integrated air compressor.

The system is designed for facility use, including within-facility transport, and should only be used under the orders of a clinician.

Prescription Use XXX  
(Part 21 CFR 801 Subpart D)

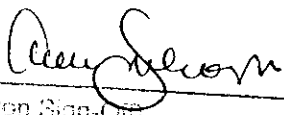
AND/OR

Over-The-Counter Use \_\_\_\_\_  
(21 CFR 807 Subpart C)

(PLEASE DO NOT WRITE BELOW THIS LINE-CONTINUE ON ANOTHER  
PAGE IF  
NEEDED)

Concurrence of CDRH, Office of Device Evaluation (ODE)

Page 1 of 1

  
Division Sign-Off  
Division of Anesthesiology, General Hospital,  
Injection Control, Dental Devices  
510(k) Number K050597