

YAHORNG

Ya Horng CO., LTD.

No. 35, Zsha Lun, Jon Zsha village,

Antin Shiang, Tainan, Taiwan, ROC

Tel: 886-6-5932201 Fax: 886-6-5935870

E-mail: lab@yahorng.com http:// www.yahorng.com

K051863
p.1/3

AUG 19 2005

"510(k) Summary"

Submitter's Name: YA HORNG Electronic Co., Ltd.

Address: *No. 35, Zaha Lun, Jon Zsha Village, Antin Shiang, 745, Taiwan, ROC*

Telephone: 886-6-5932201

FAX: 886-6-5935870

Contact Person: Dr. Jen, Ke-Min

Date Summary Prepared: July 2, 2005

Proprietary Name: Automatic Digital Wrist Blood Pressure Monitor
BP-410, BP-420

Common Name: BLOOD PRESSURE MONITOR

Classification Name: NON-INVASIVE BLOOD-PRESSURE
MEASUREMENT SYSTEM

(per 21CFR section 870.1130)

Device Class: Class II (performance standards)

Specialty: CARDIOVASCULAR

Product code: DXN

Legally Marketed *AMLUCK AUTOMATIC DIGITAL WRIST*
(Predicate) *BLOOD PRESSURE MONITOR AK-3000 /*
Device : *AK-4000*

510(k) No: K012796

YAHORNG

Ya Horng CO., LTD.

No. 35, Zsha Lun, Jon Zsha village,

Antin Shiang, Tainan, Taiwan, ROC

Tel: 886-6-5932201 Fax: 886-6-5935870

E-mail: lab@yahorng.com http:// www.yahorng.com

Description of the new device:

YA HORNG BP-410 and BP-420 use the Oscillometric method to measure the blood pressure. The Oscillometric method is adopted clinically to measure the blood pressure recently. It is not needed to use the stethoscope, as in the traditional measuring method, to monitor the Korotkov sound when deciding the systolic or diastolic pressure. The Oscillometric method senses the vibrating signal via the closed air pipe system and utilizes the microcomputer to automatically sense the characteristics of the pulse signal. Through simple calculation, the reading can reflect the accurate real blood pressure, and the systolic pressure is defined as the pressure when the cuff pressure oscillating amplitude begins to increase and the diastolic pressure as the pressure when the cuff pressure oscillating amplitude stops decreasing.

Technological Characteristics of our new device compared to the predicate device:

The technological characteristics of YA HORNG BP-410 and BP-420 are substantially equivalent to AMLUCK AK-3000 / AK-4000. There is the same Owner, AMLUNK Inc., which FDA owner number is 9040892. Especially, there are the same design specifications, the same form and intended to be used in the same manner that mean the new device BP-410 is same as the predicate device AK-4000, and the new device BP-420 is same as the predicate device AK-3000.

YAHORNG

Ya Horng CO., LTD.

No. 35, Zsha Lun, Jon Zsha village,
Antin Shiang, Tainan, Taiwan, ROC

Tel: 886-6-5932201 Fax: 886-6-5935870

E-mail: lab@yahorng.com http:// www.yahorng.com

The mainly different of the two devices are only vision appearance and all of the devices were passed the relevant EMC and Safety standards. Thus there are substantially equivalent.

Test Summary:

1. ELECTRIC SAFETY and EMC test reports,

General safety	EN 60601-1:1990+A1+A2+A11+A12+A13	PASS
	EN 1060-1:1995, EN 1060-3:1997	PASS
EMC conformity	EN 60601-1-2: 1993	PASS

2. WOVEN COTTON SHEETING

Uses the 510K Blood-Pressure Cuff

3. PERFORMANCE & CLINICAL TEST

AAMI / ANSI SP10

YA HORNG Co. Ltd. believes this information and referred document to be sufficient for the FDA to find our proposed device substantially equivalent to the predicate product and other products currently in distribution.

HSU SHENG HSIUNG

Submitter, July 2, 2005

General Manager

YA HORNG Electronic Co., Ltd.



DEPARTMENT OF HEALTH & HUMAN SERVICES

Public Health Service

Food and Drug Administration
9200 Corporate Boulevard
Rockville MD 20850

AUG 19 2005

Ya Horng Electronics Co., Ltd.
c/o Mr. Jen Ke-Min
ROC Chinese-European Industrial Research Society
No. 58, Fu-Chiun St.
Hsin-Chu City
China, (Taiwan) 300

Re: K051863

Trade Name: Automatic Digital Wrist Blood Pressure Monitor, Model BP-410 and BP-420
Regulation Number: 21 CFR 870.1130
Regulation Name: Non-Invasive Blood Pressure Measurement System
Regulatory Class: Class II (two)
Product Code: DXN
Dated: July 2, 2005
Received: July 8, 2005

Dear Mr. Ke-Min:

We have reviewed your Section 510(k) premarket notification of intent to market the device referenced above and have determined the device is substantially equivalent (for the indications for use stated in the enclosure) to legally marketed predicate devices marketed in interstate commerce prior to May 28, 1976, the enactment date of the Medical Device Amendments, or to devices that have been reclassified in accordance with the provisions of the Federal Food, Drug, and Cosmetic Act (Act) that do not require approval of a premarket approval application (PMA). You may, therefore, market the device, subject to the general controls provisions of the Act. The general controls provisions of the Act include requirements for annual registration, listing of devices, good manufacturing practice, labeling, and prohibitions against misbranding and adulteration.

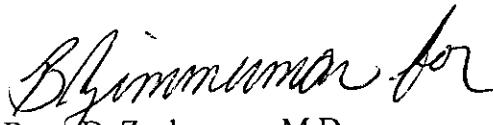
If your device is classified (see above) into either class II (Special Controls) or class III (PMA), it may be subject to such additional controls. Existing major regulations affecting your device can be found in the Code of Federal Regulations, Title 21, Parts 800 to 898. In addition, FDA may publish further announcements concerning your device in the Federal Register.

Page 2 – Mr. Jen Ke-Min

Please be advised that FDA's issuance of a substantial equivalence determination does not mean that FDA has made a determination that your device complies with other requirements of the Act or any Federal statutes and regulations administered by other Federal agencies. You must comply with all the Act's requirements, including, but not limited to: registration and listing (21 CFR Part 807); labeling (21 CFR Part 801); good manufacturing practice requirements as set forth in the quality systems (QS) regulation (21 CFR Part 820); and if applicable, the electronic product radiation control provisions (Sections 531-542 of the Act); 21 CFR 1000-1050. This letter will allow you to begin marketing your device as described in your Section 510(k) premarket notification. The FDA finding of substantial equivalence of your device to a legally marketed predicate device results in a classification for your device and thus, permits your device to proceed to the market.

If you desire specific advice for your device on our labeling regulation (21 CFR Part 801), please contact the Office of Compliance at (240) 276-0295. Also, please note the regulation entitled, "Misbranding by reference to premarket notification" (21CFR Part 807.97). You may obtain other general information on your responsibilities under the Act from the Division of Small Manufacturers, International and Consumer Assistance at its toll-free number (800) 638-2041 or (301) 443-6597 or at its Internet address <http://www.fda.gov/cdrh/industry/support/index.html>.

Sincerely yours,

A handwritten signature in black ink, appearing to read "B. Zuckerman for".

Brian D. Zuckerman, M.D.

Director

Division of Cardiovascular Devices

Office of Device Evaluation

Center for Devices and

Radiological Health

Enclosure

YAHORNG

Ya Horng CO., LTD.

No. 35, Zsha Lun, Jon Zsha village,

Antin Shiang, Tainan, Taiwan, ROC

Tel: 886-6-5932201 Fax: 886-6-5935870

E-mail: lab@yahorng.com http:// www.yahorng.com

Indications for Use

510(k) Number: K051863

Device Name: YA HORNG ELECTRONIC CO., LTD.

YA HORNG Automatic Digital Wrist Blood Pressure Monitor
BP-410, BP-420

● Indications for use:

The YA HORNG Automatic Digital Wrist Blood Pressure Monitor, Model BP-410 and BP-420, are noninvasive blood pressure measurement systems intended to measure the systolic and diastolic blood pressures and pulse rate of an adult individual, over age 18, at home by using a non-invasive technique in which an inflatable cuff is wrapped around the wrist. The cuff circumference is limited to be 5.3" – 8.5".

Prescription Use _____ AND/OR Over-The-Counter Use ✓
(Part 21 CFR 801 Subpart D) (21 CFR 807 Subpart C)

(PLEASE DO NOT WRITE BELOW THIS LINE-CONTINUE ON ANOTHER PAGE IF NEEDED)

Concurrence of CDRH, Office of Device Evaluation (ODE)

B. Himmelman
(Division Sign-Off)
Division of Cardiovascular Devices
510(k) Number K051863