

Agfa Corporation  
Premarket Notification: CR30-X

K062223

### 510(k) Summary CR30-X

Common/Classification Name: Computed Radiography, 21 CFR 892.1630

Agfa Corporation  
10 South Academy Street  
Greenville, SC 29601

SEP - 1 2006

Contact: Patrick Lynch, Prepared: July 24, 2006

#### A. LEGALLY MARKETED PREDICATE DEVICES

This is a Special 510(k) for a device modification. The modified device is Agfa's CR30-X.

The predicate device is Agfa's CR25.0 which was cleared by FDA on July 22, 2004 (K041701).

#### B. DEVICE DESCRIPTION

The predicate and newly modified devices are computed radiography imaging systems. Instead of traditional screens and photographic film for producing the diagnostic image, these systems utilize an "imaging plate," a plate coated with photo-stimulatable storage phosphors that are sensitive to X-rays and capable of retaining a latent image. After exposure, this imaging plate is inserted into a digitizer that scans it with a laser and releases the latent image in the form of light that is converted into a digital image file. The image can then be previewed on a computer workstation, adjusted if necessary then stored locally, sent to an archive, printed or sent to a softcopy capable display such as a PACS system.

The CR30-X and the CR25.0 are similar. The CR30-X utilizes an improved light collector to obtain maximum light efficiency. However, the basic principles of operation are unchanged.

#### C. INTENDED USE

The CR30-X is used to scan exposed X-ray cassettes, containing an erasable image plate (IP). This device is part of a system, consisting of X-ray cassettes with erasable phosphor image plates, an identification station for the cassettes and a workstation where the resulting digital image information is further processed and routed. It is intended that this device is only operated in a radiological environment by qualified personnel.

#### D. SUBSTANTIAL EQUIVALENCE SUMMARY

Agfa's CR30-X has the same indications for use and the same technological characteristics as the predicate device. This premarket notification has described the characteristics of the devices in sufficient detail to assure substantial equivalence. For the few characteristics that may not be precise enough to ensure equivalence, performance data was collected, and this data demonstrates substantial equivalence. In keeping with the format of a Special 510(k) for Device Modification, performance data were not included in the submission, but the

Agfa Corporation  
Premarket Notification: CR30-X

declarations provide certification that the data demonstrate equivalence.

**E. TECHNOLOGICAL CHARACTERISTICS**

The technological characteristics are the same in the proposed and predicate devices.

**F. TESTING**

The CR30-X has been tested for proper performance to specifications through various in-house reliability and imaging performance demonstration tests. The device also meets the requirements of EN 60601-1-1 and EN 60601-1-2.

**G. CONCLUSIONS**

This Special 510(k) for Device Modification submission has demonstrated Substantial Equivalence as defined and understood in the Federal Food Drug and Cosmetic Act and various guidance documents issued by the Center for Devices and Radiological Health



DEPARTMENT OF HEALTH & HUMAN SERVICES

Public Health Service

Food and Drug Administration  
10903 New Hampshire Avenue  
Document Control Room – WO66-G609  
Silver Spring, MD 20993-0002

Mr. Patrick J. Lynch  
Manager of Regulatory Affairs  
AGFA Corporation Healthcare  
10 S. Academy Street  
P.O. Box 19048  
GREENVILLE SC 29602-9048

AUG 23 2013

Re: K062223

Trade/Device Name: Agfa CR30-X  
Regulation Number: 21 CFR 892.1680  
Regulation Name: Stationary x-ray system  
Regulatory Class: II  
Product Code: MQB  
Dated: July 24, 2006  
Received: August 2, 2006

Dear Mr. Lynch:

This letter corrects our substantially equivalent letter of September 1, 2006.

We have reviewed your Section 510(k) premarket notification of intent to market the device referenced above and have determined the device is substantially equivalent (for the indications for use stated in the enclosure) to legally marketed predicate devices marketed in interstate commerce prior to May 28, 1976, the enactment date of the Medical Device Amendments, or to devices that have been reclassified in accordance with the provisions of the Federal Food, Drug, and Cosmetic Act (Act) that do not require approval of a premarket approval application (PMA). You may, therefore, market the device, subject to the general controls provisions of the Act. The general controls provisions of the Act include requirements for annual registration, listing of devices, good manufacturing practice, labeling, and prohibitions against misbranding and adulteration.

If your device is classified (see above) into class II (Special Controls), it may be subject to such additional controls. Existing major regulations affecting your device can be found in Title 21, Code of Federal Regulations (CFR), Parts 800 to 895. In addition, FDA may publish further announcements concerning your device in the Federal Register.


Please be advised that FDA's issuance of a substantial equivalence determination does not mean that FDA has made a determination that your device complies with other requirements of the Act or any Federal statutes and regulations administered by other Federal agencies. You must comply with all the Act's requirements, including, but not limited to: registration and listing (21 CFR Part 807); labeling (21 CFR Parts 801 and 809); medical device reporting (reporting of

medical device-related adverse events) (21 CFR 803); and good manufacturing practice requirements as set forth in the quality systems (QS) regulation (21 CFR Part 820). This letter will allow you to begin marketing your device as described in your Section 510(k) premarket notification. The FDA finding of substantial equivalence of your device to a legally marketed predicate device results in a classification for your device and thus, permits your device to proceed to the market.

If you desire specific advice for your device on our labeling regulation (21 CFR Parts 801 and 809), please contact the Office of *In Vitro* Diagnostic Device Evaluation and Safety at (301) 796-5450. Also, please note the regulation entitled, "Misbranding by reference to premarket notification" (21 CFR Part 807.97). For questions regarding the reporting of adverse events under the MDR regulation (21 CFR Part 803), please go to <http://www.fda.gov/MedicalDevices/Safety/ReportaProblem/default.htm> for the CDRH's Office of Surveillance and Biometrics/Division of Postmarket Surveillance.

You may obtain other general information on your responsibilities under the Act from the Division of Small Manufacturers, International and Consumer Assistance at its toll-free number (800) 638-2041 or (301) 796-7100 or at its Internet address <http://www.fda.gov/cdrh/industry/support/index.html>.

Sincerely Yours,

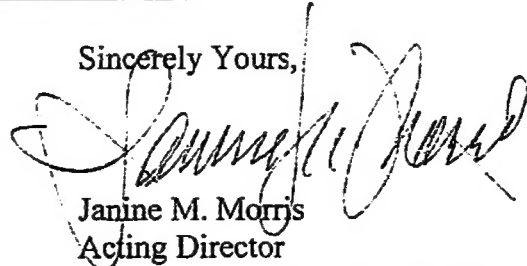A handwritten signature in black ink, appearing to read "Janine M. Morris", is written over the typed name and title.

Janine M. Morris  
Acting Director  
Division of Radiological Devices  
Office of In Vitro Diagnostic Device  
Evaluation and Safety  
Center for Devices and Radiological Health

Enclosure

Indications for Use

510(k) Number (if known):

K062223

Device Name: Agfa CR30-X

Indications for Use:

The CR30-X is indicated for use to provide diagnostic quality images to aid in physician diagnosis. The CR30-X is intended to be used mainly in chest, skeletal and gastro-intestinal x-ray imaging applications.

Prescription Use   X    
(Part 21 CFR 801 Subpart D)


AND/OR

Over-The-Counter Use             
(21 CFR 801 Subpart C)

(PLEASE DO NOT WRITE BELOW THIS LINE-CONTINUE ON ANOTHER  
PAGE IF NEEDED)

---

Concurrence of CDRH, Office of Device Evaluation (ODE)

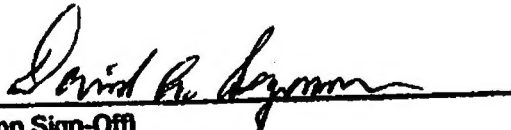

(Division Sign-Off)  
Division of Reproductive, Abdominal,  
and Radiological Devices  
510(k) Number K062223

Page   1   of   1