

**510(k) Summary**

As Required by 21 section 807.92 ( c )

- 1- Applicant Name: GENEXEL-SEIN, INC (FORMERLY SEIN ELECTRONICS)  
 2- Address: NamGyeong Building, 111, Yangjae 1-dong, Seocho-gu, Seoul, 137-891, Korea  
 3- Phone: (82) 2-575-1141  
 4- Fax: (82) 2-575-1688  
 5- Contact Person: Jea Woo Han JUN 26 2007  
 6- Date summary prepared: May 30, 2007  
 7- Official Correspondent: Mansour Consulting LLC  
 8- Address: 845 Aronson Lake Court. Roswell, GA 30075 USA  
 9- Phone: 678-908-8180  
 10- Fax: 678-623-3765  
 11- Contact Person: Jay Mansour, President  
 12- Device Trade or Proprietary Name: Full Auto Wrist Digital Blood Pressure Monitor BP1600, BP1650, BP1700 AND BP1750  
 13- Device Common or usual name: Digital Blood Pressure Monitor  
 14- Device Classification Name: Non Invasive blood pressure measuring system  
 15- Substantial Equivalency is claimed against the following device:  
 Full Auto Wrist Digital Blood Pressure Monitor, Model SE-312, manufactured by Sein Electronics Co., Ltd. (now Genexel-Sein) 510k #K012054

**16- Description of the Device:**

Full Auto Wrist Blood Pressure Monitor (Models BP1600, BP1650, BP1700 and BP1750) is intended to measure systolic and diastolic pressure and pulse rate of adults in a home care environment using wrist cuff and oscillometric method of measurement.

There are no contraindications; the subject device may be employed in the care of normotensive, hypertensive, or hypotensive patients.

The user interface panel of the device has power button, memory button and liquid crystal display ("LCD") for displaying systolic pressure, diastolic pressure and pulse rate. The device has memory capacity to store 10, 20, 60 and 90 most recent measurement results (respectively for BP1600, BP1650, BP1700 and BP1750)

The device measures blood pressure through the use of an automatic inflating cuff. Pressurization is automatically governed. If the initial inflation pressure (180mmHg) is inadequate for measurement, i.e. lower than the patient's systolic pressure, the pump will automatically re-pressurize to a preset level (30 mmHg) above the initial level. Symbols in the LCD indicate pressurization status at all times. The cuff automatically deflates with stepwise pressure drop of 5~6 mmHg per 2 pulses during blood pressure measurement.

(10)

The patient is responsible for applying the cuff, for initiating the measurements sequence by pressuring the "Power" button, and for recording results. The patient cannot alter bleed-down rate.

All system functions are preprogrammed. The user is cautioned in the instruction manual against attempting any programming or other modification.

No special training beyond basic ability to follow instruction is required. Since the products are designed for home use, detailed instructions on avoidance of practices that adversely affect the accuracy of measurements are included in the instruction manual.

**17- Intended use of the device: (refer to FDA form attached)**

To measure systolic and diastolic pressure and pulse rate of adults using wrist cuff and oscillometric method.

**18- Safety and Effectiveness of the device:**

This device is safe and effective as the predicate device cited above.  
This is better expressed in the tabulated comparison (Paragraph 14 below)

**14- Summary comparing technological characteristics with predicate device:**

Please find below a tabulated comparison supporting that this device is substantially equivalent to other medical devices in commercial distribution. Also, Equivalency overview chart path is attached. Indeed, this device is **SUBSTANTIALLY EQUIVALENT** to the predicate device. Refer to the explanations/details within the main submission.

FDA file reference number	510k # K012054
<b>TECHNOLOGICAL CHARACTERISTICS</b>	<b>Comparison result</b>
Indications for use	Identical
Target population	Identical
Design	Similar
Materials	Identical
Performance	Identical
Sterility	Not Applicable
Biocompatibility	Identical
Mechanical safety	Identical
Chemical safety	Not Applicable
Anatomical sites	Identical
Human factors	Similar
Energy used and/or delivered	Identical
Compatibility with environment and other devices	Identical
Where used	Identical
Standards met	Identical
Electrical safety	Identical

Thermal safety	Identical
Radiation safety	Not Applicable

Refer to the submission for more details.



DEPARTMENT OF HEALTH & HUMAN SERVICES

Public Health Service

Food and Drug Administration  
9200 Corporate Boulevard  
Rockville MD 20850

JUN 26 2007

GENEXEL-SEIN, INC.  
c/o Jay Mansour, MSQA, BE, LA, RAC  
Mansour Consulting, LLC  
845 Aronson Lake Court  
Roswell, Georgia 30075

Re: K071523

Trade/Device Name: Full Auto Wrist Digital Blood Pressure Monitor, Models BP1600,  
BP1650, BP1700, and BP1750

Regulation Number: 21 CFR 870.1130

Regulation Name: Non-Invasive Blood Pressure Measurement System

Regulatory Class: Class II

Product Code: DXN

Dated: May 30, 2007

Received: June 4, 2007

Dear Mr. Mansour:

We have reviewed your Section 510(k) premarket notification of intent to market the device referenced above and have determined the device is substantially equivalent (for the indications for use stated in the enclosure) to legally marketed predicate devices marketed in interstate commerce prior to May 28, 1976, the enactment date of the Medical Device Amendments, or to devices that have been reclassified in accordance with the provisions of the Federal Food, Drug, and Cosmetic Act (Act) that do not require approval of a premarket approval application (PMA). You may, therefore, market the device, subject to the general controls provisions of the Act. The general controls provisions of the Act include requirements for annual registration, listing of devices, good manufacturing practice, labeling, and prohibitions against misbranding and adulteration.

If your device is classified (see above) into either class II (Special Controls) or class III (PMA), it may be subject to such additional controls. Existing major regulations affecting your device can be found in the Code of Federal Regulations, Title 21, Parts 800 to 898. In addition, FDA may publish further announcements concerning your device in the Federal Register.

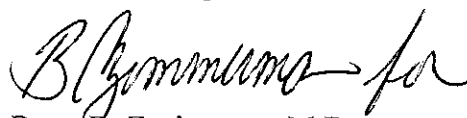
Page 2 - Mr. Jay Mansour

Please be advised that FDA's issuance of a substantial equivalence determination does not mean that FDA has made a determination that your device complies with other requirements of the Act or any Federal statutes and regulations administered by other Federal agencies. You must comply with all the Act's requirements, including, but not limited to: registration and listing (21 CFR Part 807); labeling (21 CFR Part 801); good manufacturing practice requirements as set forth in the quality systems (QS) regulation (21 CFR Part 820); and if applicable, the electronic product radiation control provisions (Sections 531-542 of the Act); 21 CFR 1000-1050.

This letter will allow you to begin marketing your device as described in your Section 510(k) premarket notification. The FDA finding of substantial equivalence of your device to a legally marketed predicate device results in a classification for your device and thus, permits your device to proceed to the market.

If you desire specific advice for your device on our labeling regulation (21 CFR Part 801), please contact the Office of Compliance at 240-276-0120. Also, please note the regulation entitled, "Misbranding by reference to premarket notification" (21CFR Part 807.97). You may obtain other general information on your responsibilities under the Act from the Division of Small Manufacturers, International and Consumer Assistance at its toll-free number (800) 638-2041 or (240) 276-3150 or at its Internet address <http://www.fda.gov/cdrh/industry/support/index.html>.

Sincerely yours,

A handwritten signature in black ink, appearing to read "B. Zuckerman for".

Bram D. Zuckerman, M.D.

Director

Division of Cardiovascular Devices

Office of Device Evaluation

Center for Devices and

Radiological Health

Enclosure

## Indications for Use

510(k) Number (if known): K071523

Device Name: Full Auto Wrist Blood Pressure Monitor, models BP1600, BP1650, BP1700 and BP1750

### Indications For Use:

To measure systolic and diastolic blood pressure and pulse rates of adults using wrist cuff and oscillometric method.

Prescription Use \_\_\_\_\_  
(Part 21 CFR 801 Subpart D)

AND/OR

Over-The-Counter Use ✓  
(21 CFR 801 Subpart C)

(PLEASE DO NOT WRITE BELOW THIS LINE-CONTINUE ON ANOTHER PAGE IF NEEDED)

---

Concurrence of CDRH, Office of In Vitro Diagnostic Devices (OIVD)

B. J. Minima  
(Division Sign-Off)

Division of Cardiovascular Devices

510(k) Number K071523

Page 1 of 1

(15)