

K071650

**ORIGINAL, TRADITIONAL 510(K) NOTIFICATION  
PERMOBIL POWERED WHEELCHAIR: C350**



Attachment 11

**510(k) Summary**

SEP 12 2007

**Submitter** Permobil AB  
Box 120  
S-861 23 Timrå  
Sweden

**Phone:** +46 60 595900

**Facsimile:** +46 60 575250

**Contact Person:** Jan Åström  
**e-mail address:** jan.astrom@permobil.se

**Date Prepared:** June, 2007

**Device name:** C350


**Classification Name:**  
Powered wheelchair

**Predicate Devices:**  
C300(Electro)( K041219) manufactured by Permobil AB.

**Intended use:**  
The intended use of the C350 powered wheelchair is to provide outdoor and indoor mobility to persons limited to a seated position that are capable of operating a powered wheelchair.

**Description of device:**  
C350 Powered Wheelchair is battery powered, rear wheel motor driven and is controlled by the PG power wheelchair VR-2 70 amp controller.  
The user interface is a joystick.  
The C350 is powered by two 12VDC 60Ah, Group M34 batteries, approximate driving range on fully charged batteries is up to 22 km (13.75 miles), depending on use and the terrain the chair is driven on.  
The chair frame is a rivet nut and welded steel construction and includes two rear drive wheels with drive units (motor, gear, brake), batteries and front pivoting casters.  
Depending on users needs, the joystick motor control is mounted to the left or right armrest.  
When the user activates the joystick, the controller receives a signal to release the brakes.  
With the brakes released, the chair is allowed to move in the direction the joystick is actuated.  
When the user releases the joystick, the chair slows to a stop and the brakes are automatically re-engaged. The solenoid electromechanical brakes allow the user stop by letting go of the joystick.

**ORIGINAL, TRADITIONAL 510(K) NOTIFICATION**  
**PERMOBIL POWERED WHEELCHAIR: C350**

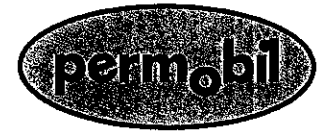

**Performance Data**

In all instances, the C350 functioned as intended.

**Substantial Equivalence**

The C350 is substantially equivalent to the C300(Electro) (#K041219). The C350 has the same intended uses and similar indications, technological characteristics and principles of operation. The minor technological differences between the Electro and its predicate device raise no new issues of safety or effectiveness. Performance data demonstrate that the C350 is as safe and effective as the C300(Electro). Thus, the C350 is substantially equivalent.



Food and Drug Administration  
9200 Corporate Boulevard  
Rockville MD 20850

Permobil AB  
% Mr. Jan Astrom  
Quality Engineer  
Box 120  
S-861 23 Timra  
Sweden

SEP 12 2007

Re: K071650  
Trade/Device Name: Powered wheelchair  
Regulation Number: 21 CFR 890.3860  
Regulation Name: Powered Wheelchair  
Regulatory Class: Class II  
Product Code: ITI  
Dated: August 9, 2007  
Received: August 17, 2007

Dear Mr. Astrom:

We have reviewed your Section 510(k) premarket notification of intent to market the device referenced above and have determined the device is substantially equivalent (for the indications for use stated in the enclosure) to legally marketed predicate devices marketed in interstate commerce prior to May 28, 1976, the enactment date of the Medical Device Amendments, or to devices that have been reclassified in accordance with the provisions of the Federal Food, Drug, and Cosmetic Act (Act) that do not require approval of a premarket approval application (PMA). You may, therefore, market the device, subject to the general controls provisions of the Act. The general controls provisions of the Act include requirements for annual registration, listing of devices, good manufacturing practice, labeling, and prohibitions against misbranding and adulteration.

If your device is classified (see above) into either class II (Special Controls) or class III (PMA), it may be subject to such additional controls. Existing major regulations affecting your device can be found in the Code of Federal Regulations, Title 21, Parts 800 to 898. In addition, FDA may publish further announcements concerning your device in the Federal Register.


Please be advised that FDA's issuance of a substantial equivalence determination does not mean that FDA has made a determination that your device complies with other requirements of the Act or any Federal statutes and regulations administered by other Federal agencies. You must comply with all the Act's requirements, including, but not limited to: registration and listing (21 CFR Part 807); labeling (21 CFR Part 801); good manufacturing practice requirements as set forth in the quality systems (QS) regulation (21 CFR Part 820); and if applicable, the electronic product radiation control provisions (Sections 531-542 of the Act); 21 CFR 1000-1050.

Page 2 – Mr. Jan Astrom

This letter will allow you to begin marketing your device as described in your Section 510(k) premarket notification. The FDA finding of substantial equivalence of your device to a legally marketed predicate device results in a classification for your device and thus, permits your device to proceed to the market.

If you desire specific advice for your device on our labeling regulation (21 CFR Part 801), please contact the Office of Compliance at (240) 276-0120. Also, please note the regulation entitled, "Misbranding by reference to premarket notification" (21CFR Part 807.97). You may obtain other general information on your responsibilities under the Act from the Division of Small Manufacturers, International and Consumer Assistance at toll-free number (800) 638-2041 or (240) 276-3150 or the Internet address <http://www.fda.gov/cdrh/industry/support/index.html>.

Sincerely yours,

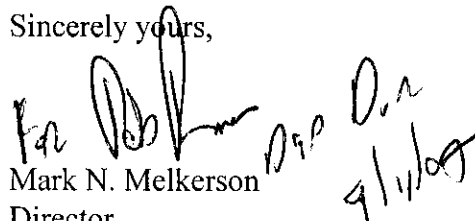A handwritten signature in black ink, appearing to read "Mark N. Melkerson". To the right of the signature, there is a date "9/11/09" and some other illegible handwritten notes.

Mark N. Melkerson  
Director  
Division of General, Restorative  
and Neurological Devices  
Office of Device Evaluation  
Center for Devices and  
Radiological Health

Enclosure

### INDICATION FOR USE

510(k) number Not assigned at the writing of this submission.

Device name: C350

**Indication for Use**

The intended use of the C350 series of the powered wheelchair is to provide indoor and outdoor mobility to persons limited to a seating position that are capable of operating a powered wheelchair.

Prescription use   X  

or


Over the counter use           

(Part 21 CFR 801 Subpart D)

(Part 21 CFR 807 Subpart C)

---

Concurrence of CDRH, Office of Device Evaluation (ODE)

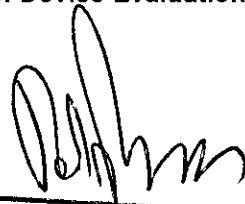

(Division Sign-Off)

Division of General, Restorative,  
and Neurological Devices

510(k) Number

14071656