

K071785

smiths

SECTION 5, 510(k) Summary

Company Information:

Smiths Medical ASD, Inc.
10 Bowman Drive
Keene, NH 03431
(603) 352-3812, prompt 4, ext 2493
Contact: Brian D. Farias
Regulatory Affairs Manager

Smiths Medical ASD, Inc.
Anesthesia and Safety Devices Division
10 Bowman Drive
Keene NH 03431 USA
Tel: +1 603 352 3812
Fax: +1 603 352 3703
www.smiths-medical.com

AUG - 2 2007

Summary Prepared: June 29, 2007

Product Name:

Trade Name: Hypodermic Needle-Pro® EDGE™ Safety Device
Hypodermic Needle-Pro® EDGE™ Safety Device with Syringe

Common Name: Hypodermic Needle with attached needle protection and Syringe with attached needle and needle protection

Classification Name: Hypodermic Single Lumen Needles (21 CFR 880.5570, Product Code FMI) and Piston Syringe (21 CFR 880.5860, Product Code FMF)

Predicate Device(s):

K041399 and K063450 (Smiths Medical ASD, Inc.) Hypodermic Needle-Pro® EDGE™ Needle Protection Device

K061194 (Smiths Medical ASD, Inc.) Portex® Needle-Pro® EDGE™ Safety Device with Syringe

K031453 (Terumo Medical Corporation) SurGuard2™ Safety Hypodermic Needle

Device Description:

This device is intended for injection or aspiration of fluids into the body. The needle protection device covers the needle after use to help prevent needle sticks. The device features a "one-piece" design of needle hub and protective sheath with a living hinge. The needle cannula is permanently affixed into the hub. The sheath has an "arrow" indicating the bevel orientation, i.e. when the sheath is oriented to the right, the bevel is in the "up position". After the procedure is completed, the needle is pressed into the sheath using a one-handed technique. As the needle enters the protective sheath, the lugs on the needle hub snap into the clips of the protective sheath resulting in a "bottom snap" and the needle is contained within the sheath. In addition to the "bottom snap", the needle protection sheath for needles equal to or less than 1" has one needle retaining hook (except 30g) and the sheath for needles longer than 1" has two hooks. The device is then discarded into a sharps container.

Indications for Use:

Hypodermic Needle-Pro® EDGE™ Safety Device: This device is intended for injection or aspiration of fluids using a Luer lock or Luer slip syringe. The needle protection device covers the needle after use to help prevent needle sticks.

Hypodermic Needle-Pro® EDGE™ Safety Device with Syringe: This device is intended for use to inject fluids into or withdraw fluids from the body. The needle protection device covers the needle after use to help prevent needle sticks.

Technological Characteristics:

The proposed and predicate devices all employ the same hinged style protective sheath that is manually activated after use.

Non-Clinical Data:

This abbreviated 510(k) submission declares conformance to the following standards:

ISO 594/1:1986(E), International Standard, *Conical fittings with a 6% (Luer) taper for syringes, needles and certain other medical equipment-Part 1: General requirements*

ISO 594-2:1998(E), International Standard, *Conical fittings with a 6% (Luer) taper for syringes, needles and certain other medical equipment-Part 2: Lock Fittings*

ISO 7864:1993(E), International Standard, *Sterile hypodermic needles for single use*

ISO 7886-1:1993, International Standard, *Sterile hypodermic syringes for single use*

Clinical Data:


Simulated clinical use studies were conducted which confirmed that the device could be used effectively, with the needle shielded inside the protection device after use.

Conclusion:

The standards compliance and simulated clinical use studies demonstrate that the proposed device is safe and effective and is substantially equivalent to the predicate devices.

Very truly yours,

SMITHS MEDICAL ASD, INC.

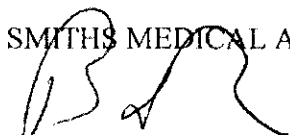

Brian D. Farias
Regulatory Affairs Manager



DEPARTMENT OF HEALTH & HUMAN SERVICES

Public Health Service

Food and Drug Administration
9200 Corporate Boulevard
Rockville MD 20850

Mr. Brian D. Farias
Regulatory Affairs Manager
Smiths Medical ASD, Incorporated
10 Bowman Drive
Keene, New Hampshire 03431

AUG - 2 2007

Re: K071785

Trade/Device Name: Hypodermic Needle-Pro® EDGE™ Safety Device

Regulation Number: 880.5570

Regulation Name: Hypodermic Single Lumen Needle

Regulatory Class: II

Product Code: FMI

Dated: June 29, 2007

Received: July 3, 2007

Dear Mr. Farias:

We have reviewed your Section 510(k) premarket notification of intent to market the device referenced above and have determined the device is substantially equivalent (for the indications for use stated in the enclosure) to legally marketed predicate devices marketed in interstate commerce prior to May 28, 1976, the enactment date of the Medical Device Amendments, or to devices that have been reclassified in accordance with the provisions of the Federal Food, Drug, and Cosmetic Act (Act) that do not require approval of a premarket approval application (PMA). You may, therefore, market the device, subject to the general controls provisions of the Act. The general controls provisions of the Act include requirements for annual registration, listing of devices, good manufacturing practice, labeling, and prohibitions against misbranding and adulteration.


If your device is classified (see above) into either class II (Special Controls) or class III (PMA), it may be subject to such additional controls. Existing major regulations affecting your device can be found in the Code of Federal Regulations, Title 21, Parts 800 to 898. In addition, FDA may publish further announcements concerning your device in the Federal Register.

Please be advised that FDA's issuance of a substantial equivalence determination does not mean that FDA has made a determination that your device complies with other requirements of the Act or any Federal statutes and regulations administered by other Federal agencies. You must comply with all the Act's requirements, including, but not limited to: registration and listing (21 CFR Part 807); labeling (21 CFR Part 801); good manufacturing practice requirements as set forth in the quality systems (QS) regulation (21 CFR Part 820); and if applicable, the electronic product radiation control provisions (Sections 531-542 of the Act); 21 CFR 1000-1050.

This letter will allow you to begin marketing your device as described in your Section 510(k) premarket notification. The FDA finding of substantial equivalence of your device to a legally marketed predicate device results in a classification for your device and thus, permits your device to proceed to the market.

If you desire specific advice for your device on our labeling regulation (21 CFR Part 801), please contact the Office of Compliance at (240) 276-0115. Also, please note the regulation entitled, "Misbranding by reference to premarket notification" (21CFR Part 807.97). You may obtain other general information on your responsibilities under the Act from the Division of Small Manufacturers, International and Consumer Assistance at its toll-free number (800) 638-2041 or (301) 443-6597 or at its Internet address <http://www.fda.gov/cdrh/industry/support/index.html>.

Sincerely yours,

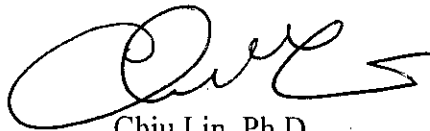A handwritten signature in black ink, appearing to read 'Chiu Lin', with a stylized flourish at the end.

Chiu Lin, Ph.D.

Director

Division of Anesthesiology, General Hospital,
Infection Control and Dental Devices

Office of Device Evaluation

Center for Devices and

Radiological Health

Enclosure

SECTION 4, Indications for Use Statement

Indications for Use

510(k) Number (if known): _____

Device Name: Hypodermic Needle-Pro® EDGE™ Safety Device

Indications for Use:

This device is intended for injection or aspiration of fluids using a Luer lock or Luer slip syringe. The needle protection device covers the needle after use to help prevent needle sticks.

Prescription Use X
(Part 21 CFR 801 Subpart D)


AND/OR

Over-The-Counter Use _____
(21 CFR 807 Subpart C)

(PLEASE DO NOT WRITE BELOW THIS LINE-CONTINUE ON ANOTHER PAGE IF NEEDED)

Concurrence of CDRH, Office of Device Evaluation (ODE)

Page 1 of 1

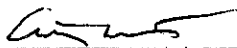

(Division Sign-Off)
Division of Anesthesiology, General Hospital,
Infection Control, Dental Devices

510(k) Number: K971785