

DEC 28 2007

K073436 page 1/2

**510(k) Summary for the  
Lutronic Corporation Spectra VRMII Laser System**

This 510(k) Summary is being submitted in accordance with the requirements of the SDMA 1990 and 21 CFR 807.92.

**1. General Information**

Submitter: Lutronic Corporation  
#403-2,3,4, Ilsan Technotown  
1141-1 Baeksok-Dong, Ilsan-Gu  
Goyang-Si, Gyeonggi-Do, 410-722  
Republic of Korea

Contact Person: Maureen O'Connell  
O'Connell Regulatory Consultants, Inc.  
5 Timber Lane  
North Reading, MA 01864  
Telephone: 978-207-1245  
Fax: 978-207-1246

Summary Preparation Date: December 20, 2007

**2. Names**

Device Name: Spectra VRMII (Q-Switched Nd:YAG)  
Laser System

Classification Name: Laser Instrument, Surgical, Powered  
Product Code: GEX  
Panel: General & Plastic Surgery

**3. Predicate Devices**

The Spectra VRMII Laser System is substantially equivalent to the Lutronic Spectra QT Q-switched Nd:YAG Laser System, the Tissue Medical Spectra-VRM Q-switched Nd:YAG Laser System, the Altus Medical CoolGlide Aesthetic Lasers, the Fotona QX Nd:YAG/KTP Laser System and the NaturaLase LP with 532 Hand Piece.

**4. Device Description**

The Spectra VRMII Laser System produces a pulsed beam of coherent near infrared (1064 nm) and visible (532nm) light. This beam is directed to the treatment zone by means of an articulated arm coupled to a handpiece. When the

beam contacts human tissue, the energy in the beam is absorbed, resulting in a very rapid, highly localized temperature increase to the target chromospheres such as melanin and tattoo particles. This increases localized temperature of the chromospheres. The instantaneous temperature increase causes fragmentation of the chromospheres to smaller particles.

By directing the beam onto specific tissue locations, using different handpieces, and controlling the treatment fluence, the intensity of the temperature of the target can be varied. The physician can optimize the effect for different applications by controlling the energy of the laser pulse and the spot size of the treatment beam.

#### **5. Indications for Use**

The Spectra VRMII Laser System is indicated for the incision, excision, ablation, vaporization of soft tissues for general dermatology, dermatologic and general surgical procedures for coagulation and hemostasis.

- 532nm Wavelength : Removal of light ink (Red, Tan, Purple, and Orange)
  - Tattoos
  - Removal of Epidermal Pigmented Lesions
  - Removal of Minor Vascular Lesions
  - Treatment of Lentigines
  - Treatment of Café-Au-Lait
  - Treatment of Seborrheic Keratoses
  - Treatment of Post Inflammatory Hyper-Pigmentation
  - Treatment of Becker's Nevi, Freckles and Nevi Spilus
- 1064nm Wavelength: Removal of dark ink (Black, Blue and Brown)
  - Tattoos
  - Removal of Nevus of Ota
  - Removal of lightening of unwanted hair with or without adjuvant preparation.
  - Treatment of Common Nevi
  - Skin resurfacing procedures for the treatment of acne scars and wrinkles

#### **6. Performance Data**

None presented.



DEPARTMENT OF HEALTH & HUMAN SERVICES

Public Health Service

Food and Drug Administration  
9200 Corporate Boulevard  
Rockville MD 20850

DEC 28 2007

Lutronic Corporation  
% O'Connell Regulatory Consultants  
Ms. Maureen O'Connell  
5 Timber Lane  
North Ridge, Massachusetts 01864

Re: K073436

Trade/Device Name: Spectra VRMII Laser System

Regulatory Number: 21 CFR 878.4810

Regulatory Name: Laser surgical instrument for use in general and plastic surgery and  
in dermatology

Regulatory Class: II

Product Code: GEX

Dated: December 5, 2007

Received: December 6, 2007

Dear Ms. O'Connell:

We have reviewed your Section 510(k) premarket notification of intent to market the device referenced above and have determined the device is substantially equivalent (for the indications for use stated in the enclosure) to legally marketed predicate devices marketed in interstate commerce prior to May 28, 1976, the enactment date of the Medical Device Amendments, or to devices that have been reclassified in accordance with the provisions of the Federal Food, Drug, and Cosmetic Act (Act) that do not require approval of a premarket approval application (PMA). You may, therefore, market the device, subject to the general controls provisions of the Act. The general controls provisions of the Act include requirements for annual registration, listing of devices, good manufacturing practice, labeling, and prohibitions against misbranding and adulteration.

If your device is classified (see above) into either class II (Special Controls) or class III (PMA), it may be subject to such additional controls. Existing major regulations affecting your device can be found in the Code of Federal Regulations, Title 21, Parts 800 to 898. In addition, FDA may publish further announcements concerning your device in the Federal Register.


Please be advised that FDA's issuance of a substantial equivalence determination does not mean that FDA has made a determination that your device complies with other requirements of the Act or any Federal statutes and regulations administered by other Federal agencies. You must comply with all the Act's requirements, including, but not limited to: registration and listing (21 CFR Part 807); labeling (21 CFR Part 801); good manufacturing practice requirements as set

forth in the quality systems (QS) regulation (21 CFR Part 820); and if applicable, the electronic product radiation control provisions (Sections 531-542 of the Act); 21 CFR 1000-1050.

This letter will allow you to begin marketing your device as described in your Section 510(k) premarket notification. The FDA finding of substantial equivalence of your device to a legally marketed predicate device results in a classification for your device and thus, permits your device to proceed to the market.

If you desire specific advice for your device on our labeling regulation (21 CFR Part 801), please contact the Center for Devices and Radiological Health’s (CDRH’s) Office of Compliance at (240) 276-0115. Also, please note the regulation entitled, “Misbranding by reference to premarket notification” (21CFR Part 807.97). For questions regarding postmarket surveillance, please contact CDRH’s Office of Surveillance and Biometric’s (OSB’s) Division of Postmarket Surveillance at (240) 276-3474. For questions regarding the reporting of device adverse events (Medical Device Reporting (MDR)), please contact the Division of Surveillance Systems at (240) 276-3464. You may obtain other general information on your responsibilities under the Act from the Division of Small Manufacturers, International and Consumer Assistance at its toll-free number (800) 638-2041 or (240) 276-3150 or at its Internet address <http://www.fda.gov/cdrh/industry/support/index.html>.

Sincerely yours,

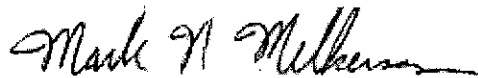

Mark N. Melkerson  
Director  
Division of General, Restorative  
and Neurological Devices  
Office of Device Evaluation  
Center for Devices and  
Radiological Health

Enclosure

**Indications for Use**

510(k) Number (if known): K073436

Device Name: Spectra VRMII Laser System

**Indications for Use:**

The Spectra VRMII Laser System is indicated for the incision, excision, ablation, vaporization of soft tissues for general dermatology, dermatologic and general surgical procedures for coagulation and hemostasis.

- 532nm Wavelength : Removal of light ink (Red, Tan, Purple, and Orange)
  - Tattoos
  - Removal of Epidermal Pigmented Lesions
  - Removal of Minor Vascular Lesions
  - Treatment of Lentigines
  - Treatment of Café-Au-Lait
  - Treatment of Seborrheic Keratoses
  - Treatment of Post Inflammatory Hyper-Pigmentation
  - Treatment of Becker's Nevi, Freckles and Nevi Spilus
- 1064nm Wavelength: Removal of dark ink (Black, Blue and Brown)
  - Tattoos
  - Removal of Nevus of Ota
  - Removal of lightening of unwanted hair with or without adjuvant preparation.
  - Treatment of Common Nevi
  - Skin resurfacing procedures for the treatment of acne scars and wrinkles

Prescription Use X  
(Part 21 CFR 801 Subpart D)

AND/OR

Over The Counter Use \_\_\_\_\_  
(Part 21 CFR 801 Subpart C)

(PLEASE DO NOT WRITE BELOW THIS LINE – CONTINUE ON ANOTHER PAGE IF NEEDED)

Concurrence of CDRH, Office of Device Evaluation (ODE)

  
(Division Sign-Off)

Division of General, Restorative,  
and Neurological Devices

510(k) Number K073436