

## 510(k) SUMMARY

MAR - 5 2008

This summary of 510(k) safety and effectiveness information is being submitted in accordance with the requirements of The Safe Medical Devices Act of 1990 (SMDA 1990) and 21 CFR 807.92.

The assigned 510(k) number is: K080322

1. *Name of Submitter, Contact Person and Date Summary Prepared:*

P. Narayan Nayak  
Director, Systems Development  
Hycor Biomedical Inc.  
7272 Chapman Avenue  
Garden Grove, CA 92841  
Phone: (714) 933-3145  
Fax: (714) 891-4153

Summary Prepared On: February 1, 2008

2. *Device Name:*

Trade/Proprietary Name: The HY•TEC™ Extended Specific IgE EIA, MCS  
Assay Using The Tecan Freedom EVO® RSP 200

Common/Usual Name: Enzyme Immunoassay System for Specific IgE

Classification Name: Radioallergosorbent (RAST) Immunological Test  
System

3. *Legally Marketed Equivalent Device Name:*

We are claiming substantial equivalence to the The HY•TEC™ Automated EIA System for Total IgE and Specific IgE, MCS Assay Using The HY•TEC™ 480, cleared 510(k) K941278.

4. *Intended Use of the Device:*

The HY•TEC™ Extended Specific IgE EIA, MCS Assay Using The Tecan Freedom EVO® RSP 200 is an enzyme immunoassay (EIA) method for the quantitative determination of allergen-specific IgE concentrations in human serum. The assay system is for in-vitro diagnostic use.

5. *Description of the Device:*

The EVO RSP 200 Automated EIA instrument is a fully-automated system, which performs sample dilution and pipetting, incubation, washing, reading and data analysis and prints reports. The EVO RSP 200 Automated EIA System for Specific IgE performs the same functions as the HY•TEC 480 Automated EIA System. The specific IgE assay for both systems, the HY•TEC 480 and the EVO RSP 200, is the same (utilizes same reagents and procedures) except for the following:

A zero calibrator is being added to the current set of five calibrators (0.35, 0.7, 3.5, 17.5 and 100 IU/mL) for the EVO Specific IgE assay.

The incubation times and temperatures have been changed.

The HY•TEC 480 and EVO RSP 200 allergy systems are standardized using anti-IgE discs and reference sera calibrated against WHO 2<sup>nd</sup> IRP 75/502 and offers a broad menu of specific allergens. An allergen-coated paper disc is incubated with a serum sample. Non-specific proteins are removed by washing and the disc is incubated with enzyme-labeled monoclonal anti human IgE conjugate. Following a second wash, substrate color is developed. The results are read spectrophotometrically against a calibration curve; results are reported in both IU/mL and Classes.

6. *Device Comparison, Verification and Results:*

Validation of software, instrument functions and assay performance, including correlation with the predicate, intra and inter assay precision, limits of detection and quantitation, and dilution linearity demonstrate the acceptability of the EVO RSP 200.

7. *Conclusion:*

The EVO RSP 200 has the same technological characteristics and intended use as the HY•TEC 480. Therefore, the modified device does not raise any new safety or effectiveness issues.



DEPARTMENT OF HEALTH & HUMAN SERVICES

Public Health Service

Food and Drug Administration  
2098 Gaither Road  
Rockville MD 20850

Hycor Biomedical, Inc.  
c/o Mr. P. Narayan Nayak  
Director, Systems Development  
7272 Chapman Ave.  
Garden Grove, CA 92841

MAR - 5 2008

Re: k080322

Trade/Device Name: HY•TEC™ Extended Specific IgE EIA, MCS Assay using Tecan  
Freedom EVO® RSP 200

Regulation Number: 21 CFR 866.5750

Regulation Name: Radioallergosorbent (RAST) immunological test system

Regulatory Class: Class II

Product Code: DHB

Dated: February 4, 2008

Received: February 6, 2008

Dear Mr. Nayak:

We have reviewed your Section 510(k) premarket notification of intent to market the device referenced above and have determined the device is substantially equivalent (for the indications for use stated in the enclosure) to legally marketed predicate devices marketed in interstate commerce prior to May 28, 1976, the enactment date of the Medical Device Amendments, or to devices that have been reclassified in accordance with the provisions of the Federal Food, Drug, and Cosmetic Act (Act) that do not require approval of a premarket approval application (PMA). You may, therefore, market the device, subject to the general controls provisions of the Act. The general controls provisions of the Act include requirements for annual registration, listing of devices, good manufacturing practice, labeling, and prohibitions against misbranding and adulteration.

If your device is classified (see above) into either class II (Special Controls) or class III (PMA), it may be subject to such additional controls. Existing major regulations affecting your device can be found in the Code of Federal Regulations, Title 21, Parts 800 to 898. In addition, FDA may publish further announcements concerning your device in the Federal Register.

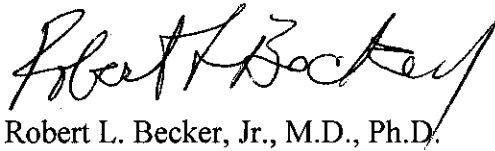
Please be advised that FDA's issuance of a substantial equivalence determination does not mean that FDA has made a determination that your device complies with other requirements of the Act or any Federal statutes and regulations administered by other Federal agencies. You must comply with all the Act's requirements, including, but not limited to: registration and listing (21 CFR Part 807); labeling (21 CFR Part 801); good manufacturing practice requirements as set forth in the quality systems (QS) regulation (21 CFR Part 820). This letter will allow you to begin marketing your device as described in your Section 510(k) premarket notification. The

Page 2 –

FDA finding of substantial equivalence of your device to a legally marketed predicate device results in a classification for your device and thus, permits your device to proceed to the market.

If you desire specific information about the application of labeling requirements to your device, or questions on the promotion and advertising of your device, please contact the Office of In Vitro Diagnostic Device Evaluation and Safety at (240) 276-0450. Also, please note the regulation entitled, "Misbranding by reference to premarket notification" (21 CFR Part 807.97). For questions regarding postmarket surveillance, please contact CDRH's Office of Surveillance and Biometric's (OSB's) Division of Postmarket Surveillance at 240-276-3474. For questions regarding the reporting of device adverse events (Medical Device Reporting (MDR)), please contact the Division of Surveillance Systems at 240-276-3464. You may obtain other general information on your responsibilities under the Act from the Division of Small Manufacturers, International and Consumer Assistance at its toll-free number (800) 638-2041 or (240) 276-3150 or at its Internet address <http://www.fda.gov/cdrh/industry/support/index.html>.

Sincerely yours,

A handwritten signature in black ink, appearing to read "Robert L. Becker, Jr.", with a stylized flourish at the end.

Robert L. Becker, Jr., M.D., Ph.D.

Director

Division of Immunology and Hematology Devices

Office of In Vitro Diagnostic Device Evaluation and Safety

Center for Devices and Radiological Health

Enclosure

## 5.0 Indications for Use Statement

### INDICATIONS FOR USE STATEMENT

510(k) Number (if known): K080322

Device Name: The HY•TEC™ Extended Specific IgE EIA, MCS Assay Using  
The Tecan Freedom EVO® RSP 200

#### Indications For Use:

The HY•TEC™ Extended Specific IgE EIA, MCS Assay Using The Tecan Freedom EVO® RSP 200 is an enzyme immunoassay (EIA) method for the quantitative determination of allergen-specific IgE concentrations in human serum. The assay system is for in-vitro diagnostic use. Measurement of specific allergen antibodies may aid in the diagnosis of asthma, allergies, and other pulmonary disorders.

Prescription Use X  
(Part 21 CFR 801 Subpart D)

OR

Over-The-Counter Use \_\_\_\_\_  
(21 CFR 801 Subpart C)

(PLEASE DO NOT WRITE BELOW THIS LINE-CONTINUE ON ANOTHER PAGE IF  
NEEDED)

Concurrence of CDRH, Office of In Vitro Diagnostic Device Evaluation and Safety (OIVD)

Maria M Chen  
Division Sign-Off  
Office of In Vitro Diagnostic Device  
Evaluation and Safety

510(k) K080322