

10091440

510(k) SUMMARY (per CFR21 807.92(c))

GENERAL INFORMATION:

OCT - 9 2009

**510k Owner's Name
Address**

Bovie Medical
7100 30th Avenue North
St. Petersburg, Florida 33710-2902

Contact Person

Richard A. Kozloff
Vice-President; Quality Assurance/Regulatory Affairs
Telephone #: (727) 384-2323
FAX Number: (727) 347-9144

Date Prepared:

May 12, 2009

DEVICE DESCRIPTION:

Trade Name:

Laparoscopic Saline Enhanced Electrosurgical Resection
(SEER) Device

Common Name:

Resection (Cutting and Coagulation) Device

Classification Name:

Electrosurgical Cutting and Coagulation Devices and
Accessories (21CFR 878.4400; Class II;
Product Code: GEI)

Predicate Device:

Bovie Medical Resection Device (K-082568)

510(k) SUMMARY (per CFR21 807.92(c))

INTENDED USE:

The Resection Device is a sterile, single use electrosurgical device intended to be used in conjunction with an electrosurgical generator for the laparoscopic delivery of radiofrequency ("RF") current and sterile saline for cutting and coagulating soft tissue.

DEVICE COMPONENTS AND OPERATION:

1) Insulated Handle:

The Insulated Handle encases the controlling mechanism for the flow of saline, and activation of the RF current for the device.

The activation of RF current is accomplished by a single push button on the top of the handle.

The Handle has a flow control mechanism so the flow of saline can be regulated by the user within the sterile field. The tubing length is approximately ten (10) feet in length and incorporates an I.V. spike on the end to attach directly to a hanging IV (saline) bag.

The Handle power cord is approximately ten (10) feet in length and incorporates a 3-prong electrical plug.

The insulation of the Insulated Handle and Power Cord meets the requirements for Dielectric Withstands of Accessories.

2) Shaft and Electrode Tip:

The electrode tip delivers RF energy for cutting and coagulation and delivers saline which is gravity-fed from an Intravenous bag to the tip.

The electrode shaft measures 32 centimeters (12.6 inches) in length and is intended to be placed through a cannula with a minimum diameter of 10 millimeters.

510(k) SUMMARY (per CFR21 807.92(c))

This device use technology substantially equivalent to the Bovie Medical Resection Device (K-082568). Both consist of a sintered stainless steel electrode tip that is used to cut and coagulate tissue through the utilization of high frequency radiofrequency energy. The difference between the two devices is that the predicate device is used for open surgical procedures and the Laparoscopic Resection Device is used for laparoscopic procedures.

Laparoscopic Resection Devices are provided sterile, sterilized using ethylene oxide gas, and are for single use only.

The Laparoscopic Resection Device conforms to the requirements of safety standard IEC 60601-2-2.

There are no significant differences in technology, intended use, or performance between the Laparoscopic Resection Device and the given predicate device. There are no new questions raised regarding safety or effectiveness.



DEPARTMENT OF HEALTH & HUMAN SERVICES

Food and Drug Administration
10903 New Hampshire Avenue
Document Mail Center - WO66-G609
Silver Spring, MD 20993-0002

Bovie Medica Inc.
% Mr. Richard Kozloff
7100 30th Avenue North
Saint Petersburg, Florida 33710

OCT - 9 2009

Re: K091440

Trade/Device Name: Laparoscopic Resection Device, model sr326L
Regulation Number: 21 CFR 878.4400
Regulation Name: Electrosurgical cutting and coagulation device and accessories
Regulatory Class: II
Product Code: GEI
Dated: August 21, 2009
Received: August 25, 2009

Dear Mr. Kozloff:

We have reviewed your Section 510(k) premarket notification of intent to market the device referenced above and have determined the device is substantially equivalent (for the indications for use stated in the enclosure) to legally marketed predicate devices marketed in interstate commerce prior to May 28, 1976, the enactment date of the Medical Device Amendments, or to devices that have been reclassified in accordance with the provisions of the Federal Food, Drug, and Cosmetic Act (Act) that do not require approval of a premarket approval application (PMA). You may, therefore, market the device, subject to the general controls provisions of the Act. The general controls provisions of the Act include requirements for annual registration, listing of devices, good manufacturing practice, labeling, and prohibitions against misbranding and adulteration.

If your device is classified (see above) into either class II (Special Controls) or class III (PMA), it may be subject to additional controls. Existing major regulations affecting your device can be found in the Code of Federal Regulations, Title 21, Parts 800 to 898. In addition, FDA may publish further announcements concerning your device in the Federal Register.

Please be advised that FDA's issuance of a substantial equivalence determination does not mean that FDA has made a determination that your device complies with other requirements of the Act or any Federal statutes and regulations administered by other Federal agencies. You must comply with all the Act's requirements, including, but not limited to: registration and listing (21 CFR Part 807); labeling (21 CFR Part 801); medical device reporting (reporting of medical device-related adverse events) (21 CFR 803); good manufacturing practice requirements as set forth in the quality systems (QS) regulation (21 CFR Part 820); and if applicable, the electronic product radiation control provisions (Sections 531-542 of the Act); 21 CFR 1000-1050.

Page 2- Mr. Richard Kozloff

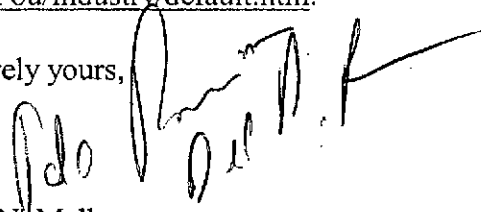
If you desire specific advice for your device on our labeling regulation (21 CFR Part 801), please go to <http://www.fda.gov/AboutFDA/CentersOffices/CDRH/CDRHOffices/ucm115809.htm> for the Center for Devices and Radiological Health's (CDRH's) Office of Compliance. Also, please note the regulation entitled, "Misbranding by reference to premarket notification" (21CFR Part 807.97). For questions regarding the reporting of adverse events under the MDR regulation (21 CFR Part 803), please go to

<http://www.fda.gov/MedicalDevices/Safety/ReportaProblem/default.htm> for the CDRH's Office of Surveillance and Biometrics/Division of Postmarket Surveillance.

You may obtain other general information on your responsibilities under the Act from the Division of Small Manufacturers, International and Consumer Assistance at its toll-free number (800) 638-2041 or (301) 796-7100 or at its Internet address

<http://www.fda.gov/MedicalDevices/ResourcesforYou/Industry/default.htm>.

Sincerely yours,

A handwritten signature in black ink, appearing to read 'Mark N. Melkerson', is written over the typed name and title.

Mark N. Melkerson

Director

Division of Surgical, Orthopedic,
and Restorative Devices

Office of Device Evaluation

Center for Devices and

Radiological Health

Enclosure

Indications for Use

510(k) Number (if known): K091440

Device Name: Laparoscopic Resection Device

Indications for Use:

The Laparoscopic Resection Device is a sterile, single use electrosurgical device intended to be used in conjunction with an electrosurgical generator for the laparoscopic delivery of radiofrequency ("RF") current and sterile saline for cutting and coagulating soft tissue. The device is not intended for contraceptive tubal coagulation (permanent female sterilization).

Prescription Use ✓
(Part 21 CFR 801 Subpart D)

AND/OR

Over-The-Counter Use
(21 CFR 801 Subpart C)

(PLEASE DO NOT WRITE BELOW THIS LINE - CONTINUE ON ANOTHER PAGE
IF NEEDED)

Concurrence of CDRH, Office of Device Evaluation (ODE)



(Division Sign-Off)

Division of Surgical, Orthopedic,
and Restorative Devices

510(k) Number K091440